

SMARTX: Porewater Sulfide

July 2026 Samples

2026-07-24

Code Set up

run information

```
##Before or After Reruns (YES or NO)
After_Reruns = "NO"
#If NO will not write final data file
#And Will Write a Rerun File

###things that need to be changed
sample_month = "July"
sample_year = "2026"
Run_Date = as.Date("07/23/2026", format = "%m/%d/%Y")
Run_by = "Melanie Giessner" #Instrument user
Script_run_by = "Melanie Giessner" #Code user

Run_notes="440-80 & 450-80 not collected" #any notes from run

#Import Plate Data

folder_path_rawdata <- paste0("Raw Data")
Raw_plates <- list.files(path = folder_path_rawdata, full.names = TRUE, pattern = ("(202607|202608)(.+Plate)"))
print(Raw_plates)

## [1] "Raw Data/20260720_SMARTX_H2S_Plate1.csv"
## [2] "Raw Data/20260720_SMARTX_H2S_Plate2.csv"
## [3] "Raw Data/20260720_SMARTX_H2S_Plate3.csv"
## [4] "Raw Data/20260720_SMARTX_H2S_Plate4.csv"
## [5] "Raw Data/20260720_SMARTX_H2S_Plate5.csv"
## [6] "Raw Data/20260722_SMARTX_H2S_Plate6.csv"
## [7] "Raw Data/20260722_SMARTX_H2S_Plate7.csv"
## [8] "Raw Data/20260722_SMARTX_H2S_Plate8.csv"
## [9] "Raw Data/20260722_SMARTX_H2S_Plate9.csv"

# Define the file path for QAQC log file - NO Need to change just check year
log_path <- here("Porewater QAQC Logs", "GCRw_H2S_QAQC_Log.csv")

# Define final path for data be sure to change year and month and run number
final_path <- "Processed Data/GCRw_SMARTX_Porewater_H2S_202605.csv"
rerun_path <- "Processed Data/GCRw_SMARTX_Porewater_H2S_202605_reruns.csv"



##Read in metadata and create similar sample IDs for matching to samples
##Read in data
```


```
##Remove High CV Stds
```

```
##Check R2 & Remove Standards if Necessary
```

```
## Joining with 'by = join_by(Plate, IDs)'
```

Plot Standard Curves

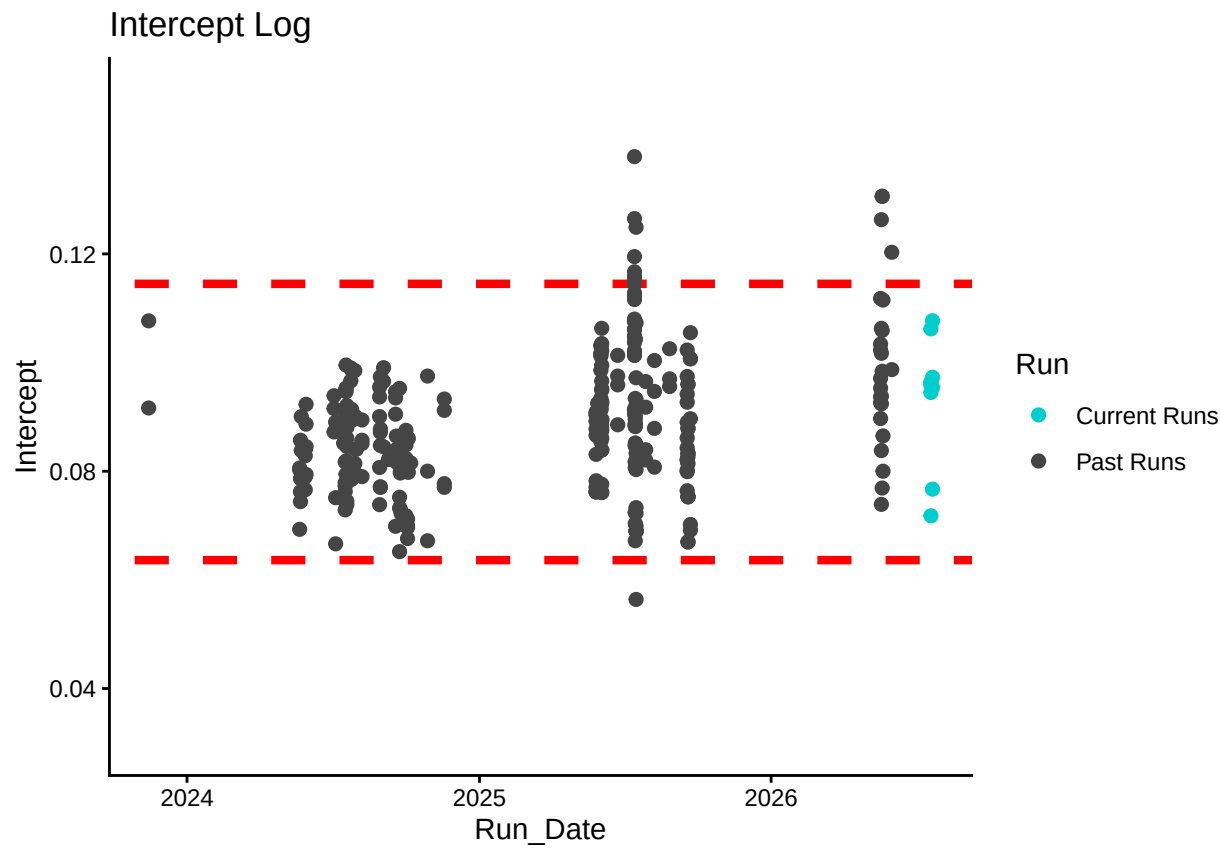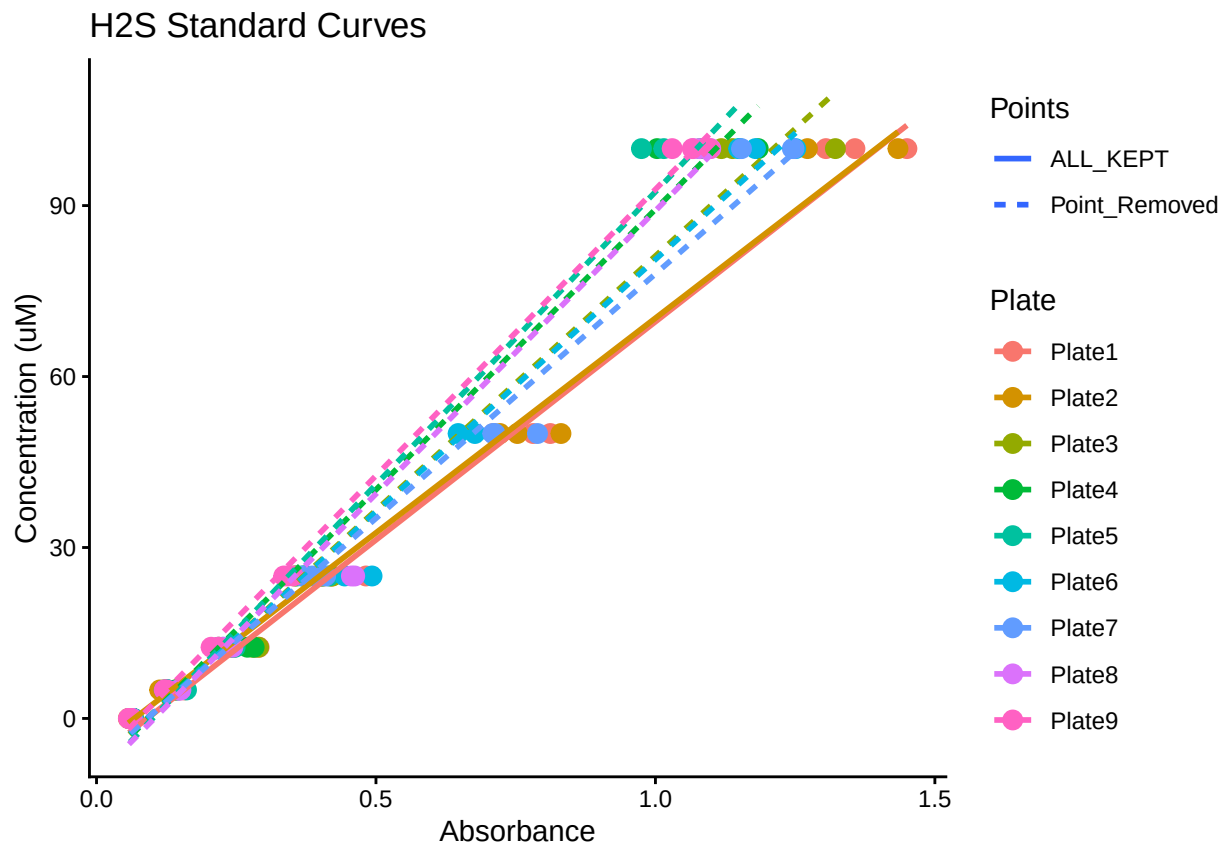

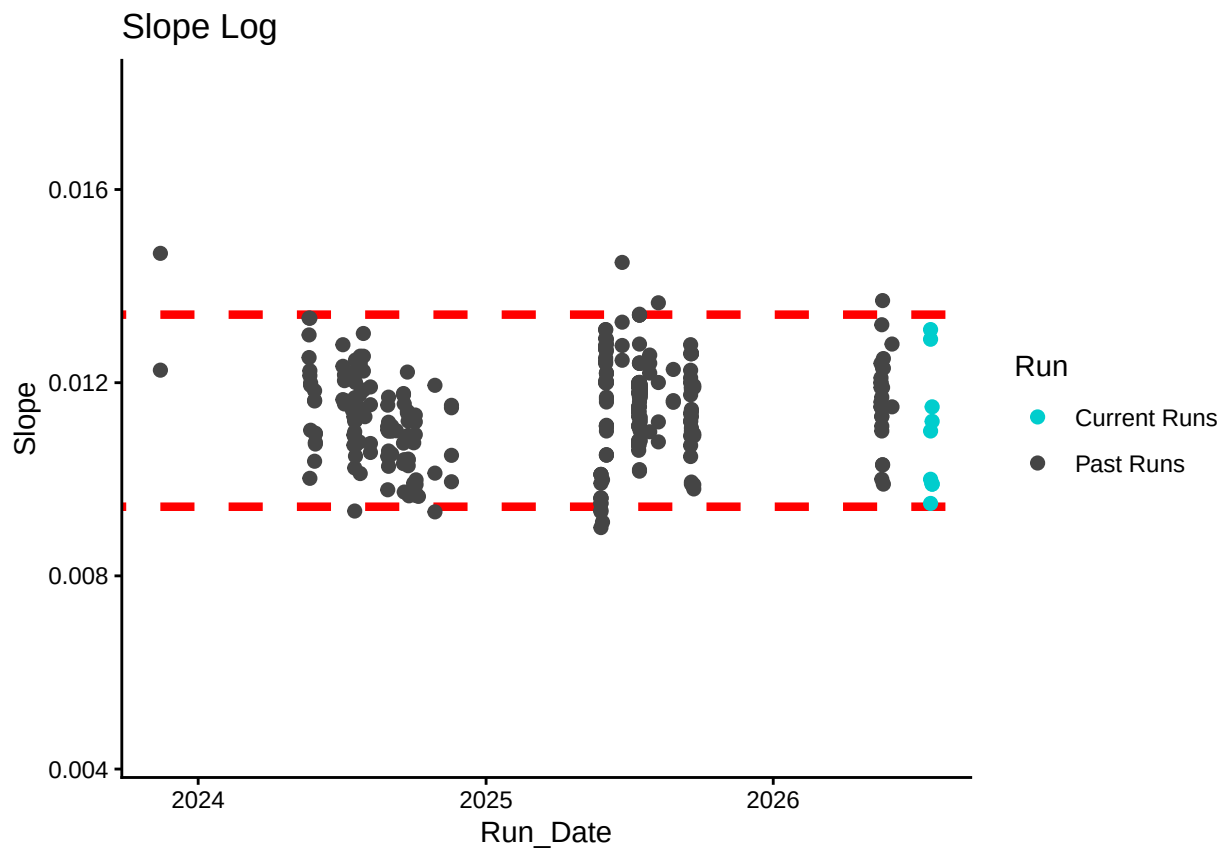


Return Flags For Standard Curve


Plate	R2_Flag	CV_Flag
Plate1	R2_Good	CV_Good
Plate2	R2_Good	CV_Good
Plate3	R2_Bad	CV_Good
Plate3	R2_Bad	CV_High
Plate4	R2_Bad	CV_Good
Plate5	R2_Bad	CV_Good
Plate5	R2_Bad	CV_High
Plate6	R2_Good	CV_Good
Plate7	R2_Good	CV_Good
Plate8	R2_Bad	CV_Good
Plate9	R2_Good	CV_Good

Plate	Curve_RMVD
Plate1	KEPT
Plate2	KEPT
Plate3	REMOVED
Plate4	REMOVED
Plate5	REMOVED
Plate6	KEPT
Plate7	KEPT
Plate8	REMOVED
Plate9	KEPT

Method Detection Limit

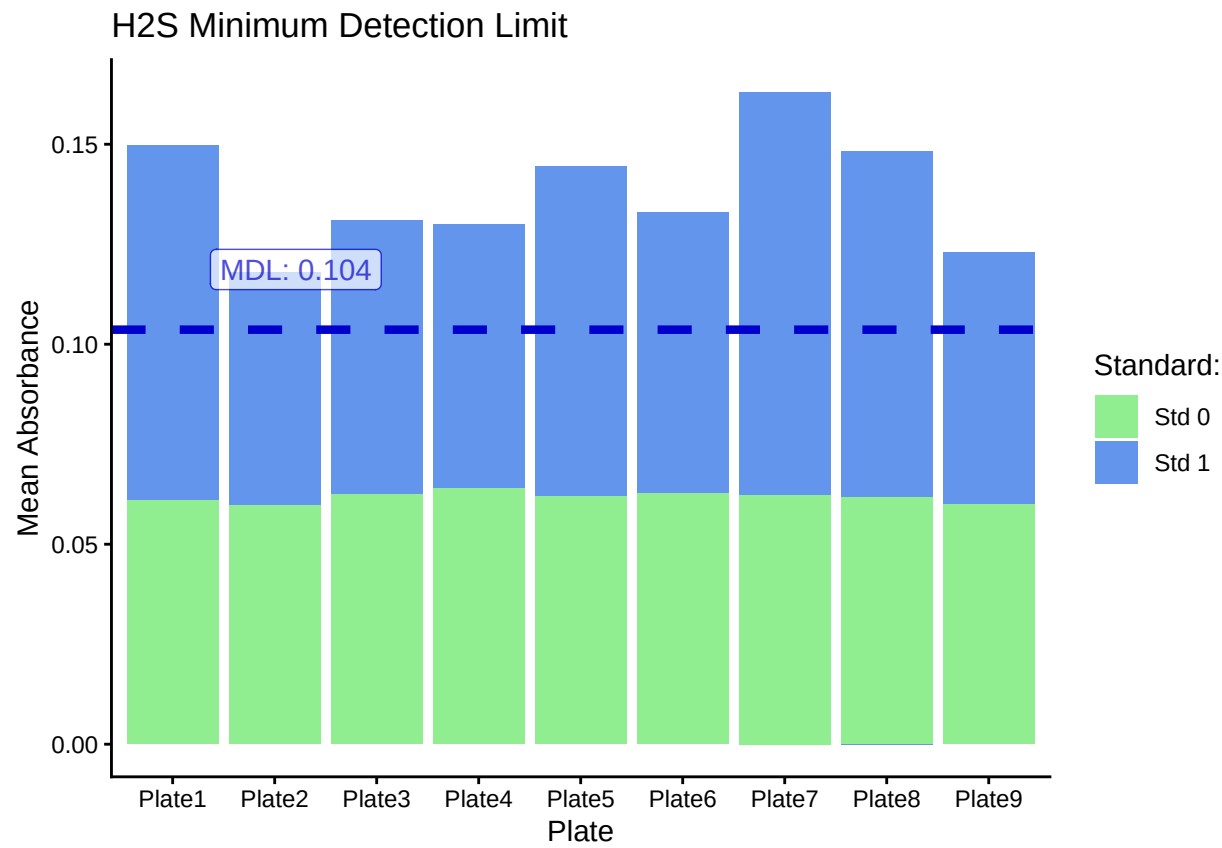

Compare Check Standards to Standards

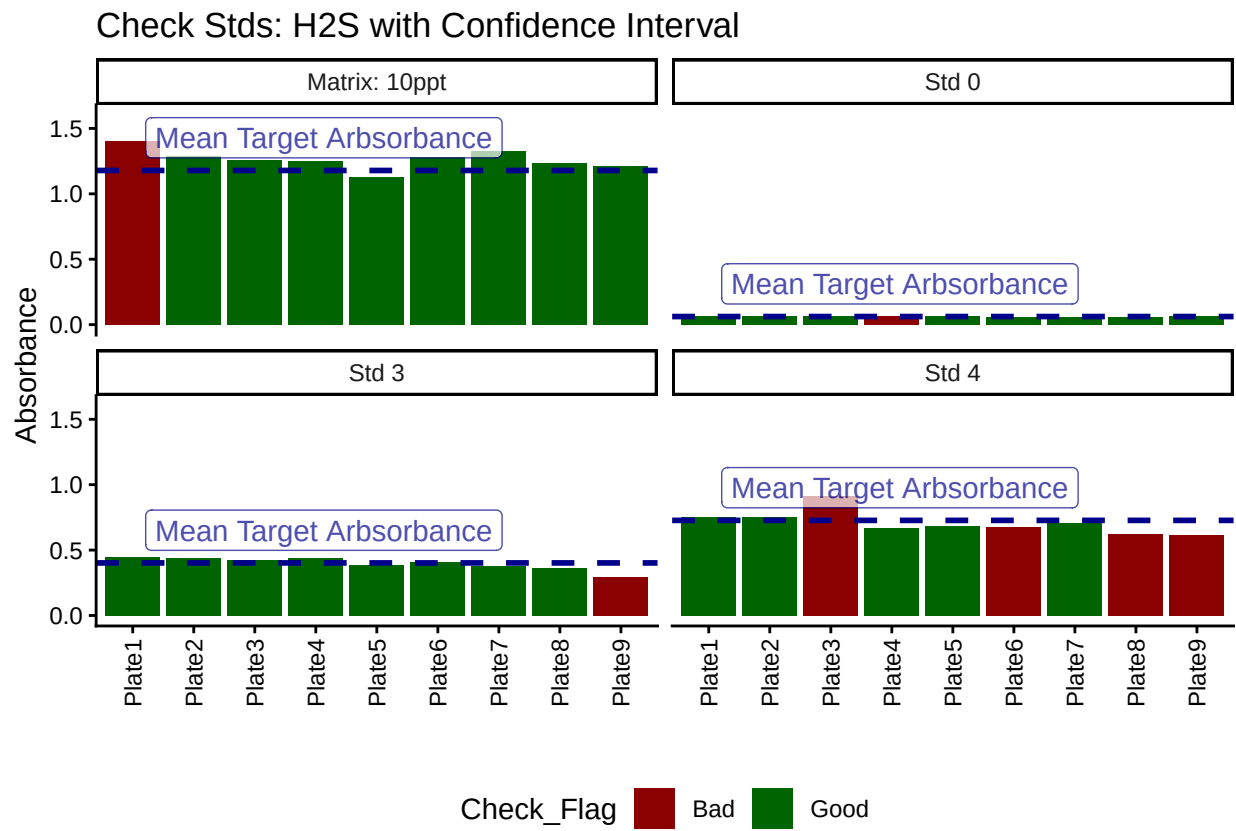

Display Any Check Flags


Plate	CHKs_Flag
Plate1	Std5Chk_Bad
Plate3	Std4Chk_Bad
Plate4	Std0Chk_Bad
Plate6	Std4Chk_Bad
Plate8	Std4Chk_Bad
Plate9	Std3Chk_Bad
Plate9	Std4Chk_Bad

Plate	CHK_CV_Flags
Plate1	Std4Chk_CV_High
Plate3	Std4Chk_CV_High


##Calculate Sulfide Concentrations & Add Concentration Flags
##Remove High CV Samples and Average

Plot Samples


Sample triplicate means and sd dev before bad reps removed

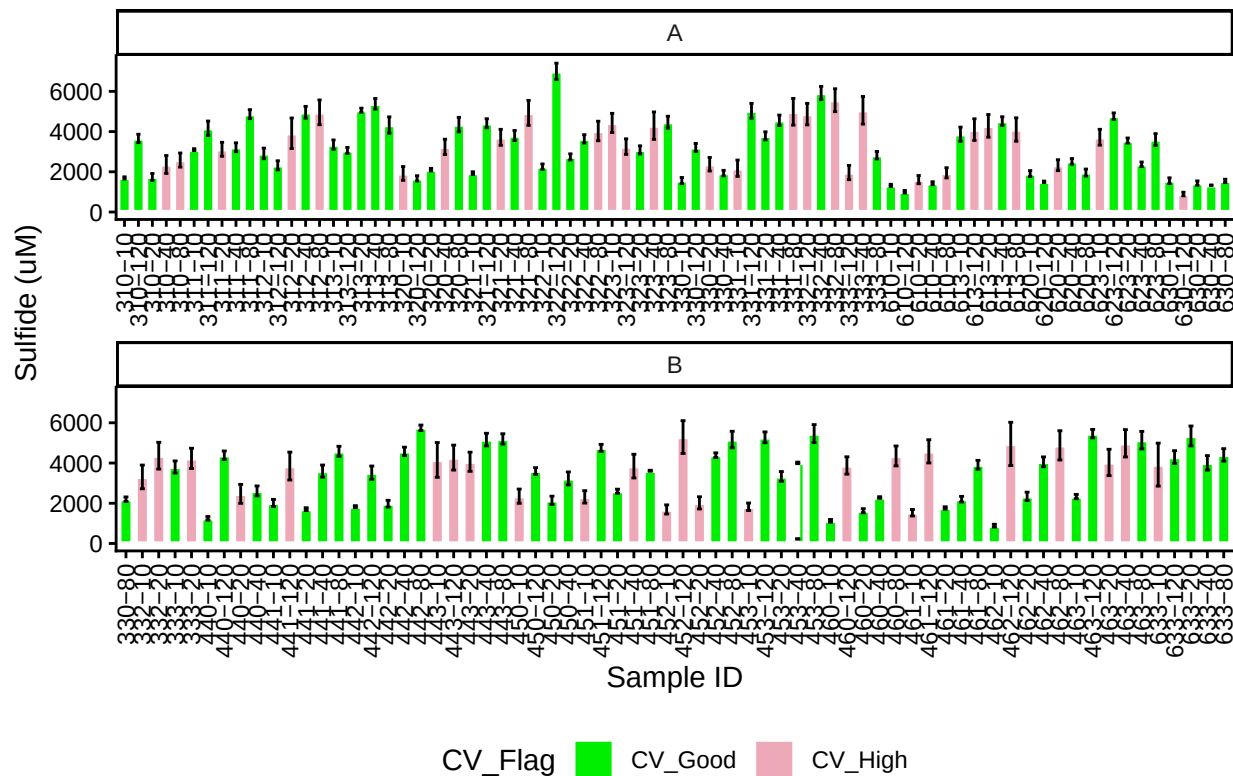

Sample triplicate means and sd dev after bad reps removed

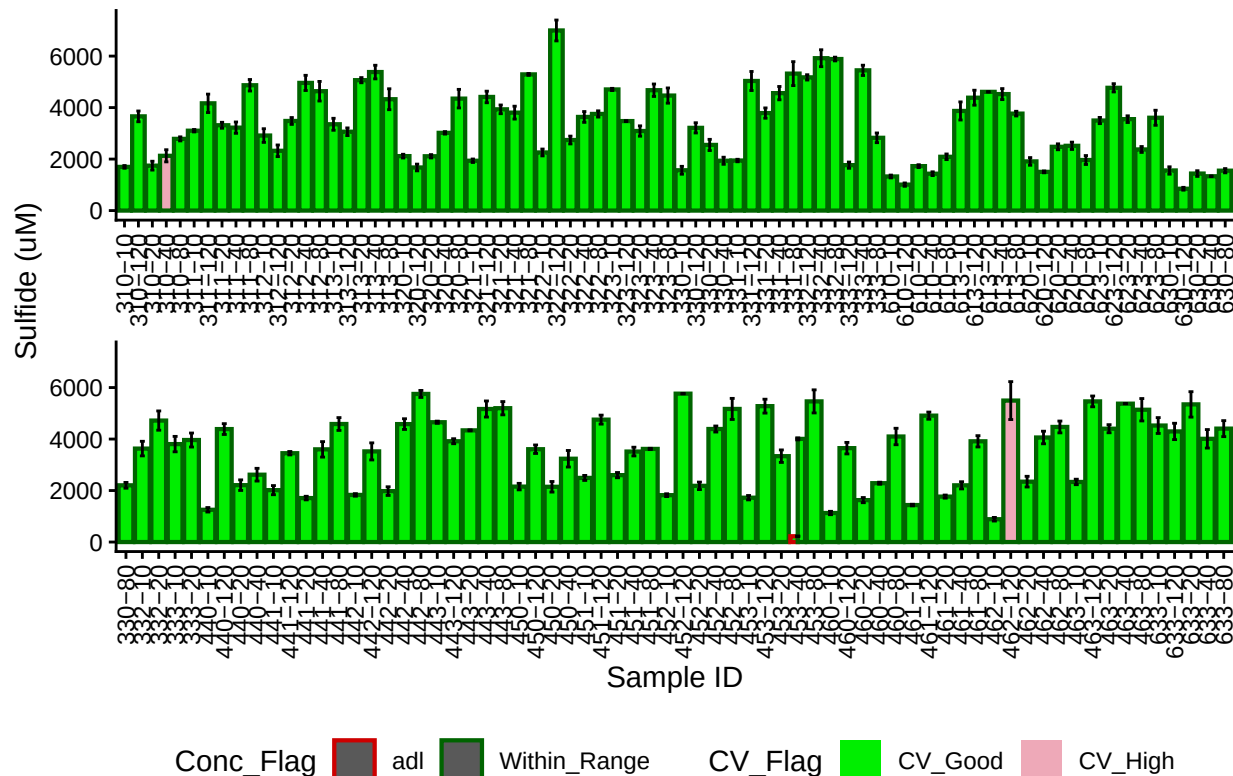

Duplicate QAQC

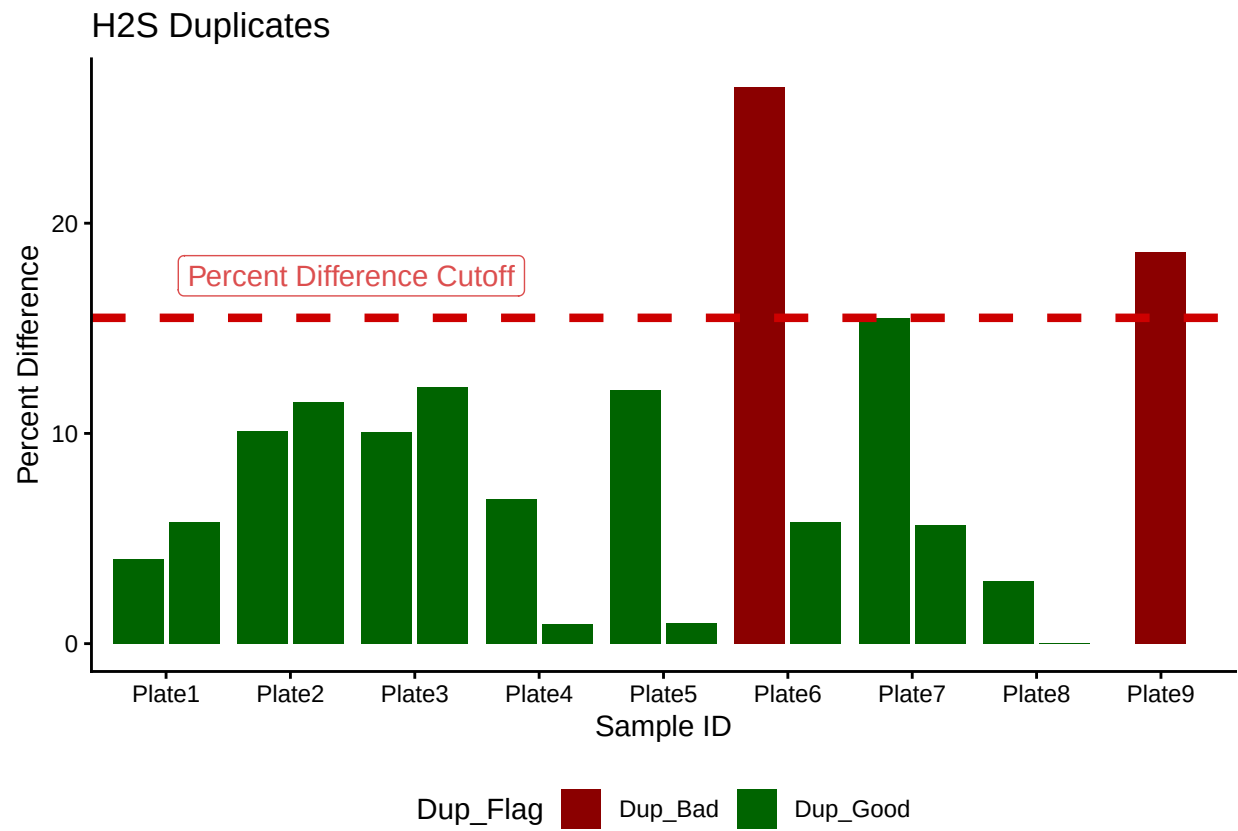


Display Dup Flags

[1] ">60% Dups Pass"

Plate	Dup_Flag
Plate1	Dup_Good
Plate2	Dup_Good
Plate3	Dup_Good
Plate4	Dup_Good
Plate5	Dup_Good
Plate6	Dup_Bad
Plate7	Dup_Good
Plate8	Dup_Good
Plate9	Dup_Bad

Spike QAQC

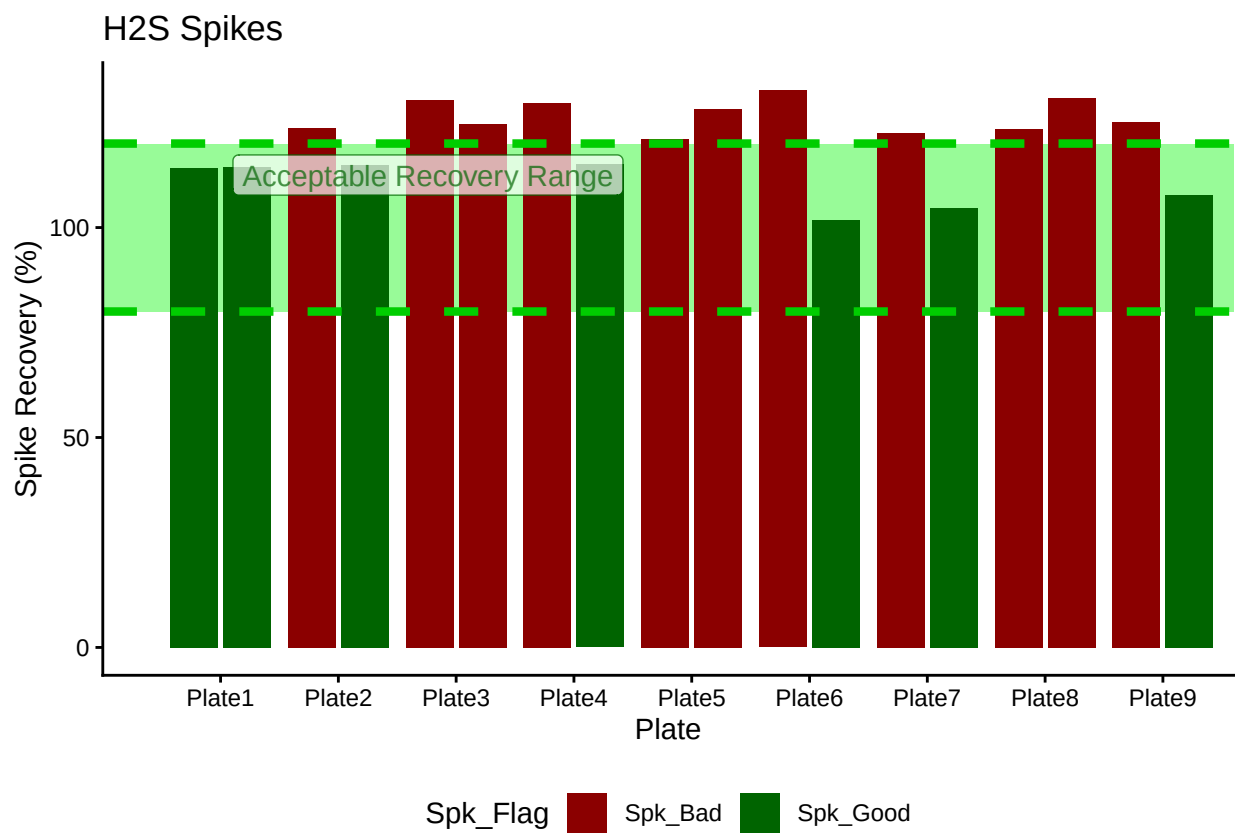

Display Spike Flags

```
## [1] "Spikes Paired"
```

Plate	Spike_Flag
Plate1	Spk_Good
Plate2	Spk_Bad
Plate3	Spk_Bad
Plate4	Spk_Bad
Plate5	Spk_Bad
Plate6	Spk_Bad
Plate7	Spk_Bad
Plate8	Spk_Bad
Plate9	Spk_Bad

```
## [1] "<60% of Spikes Pass, Rerun"
```

Merge samples with Metadata & check all samples are present

```
## Some sample IDs are missing from data.
```

```
## [1] "440-80" "450-80"
```

```
##Add Flags
```


```
##Organize Data
```

```
## Adding missing grouping variables: 'Sample_ID.x'  
## Adding missing grouping variables: 'Sample_ID.x'
```


Visualize Data

```
## Visualize Data
```

C3: Porewater H2S

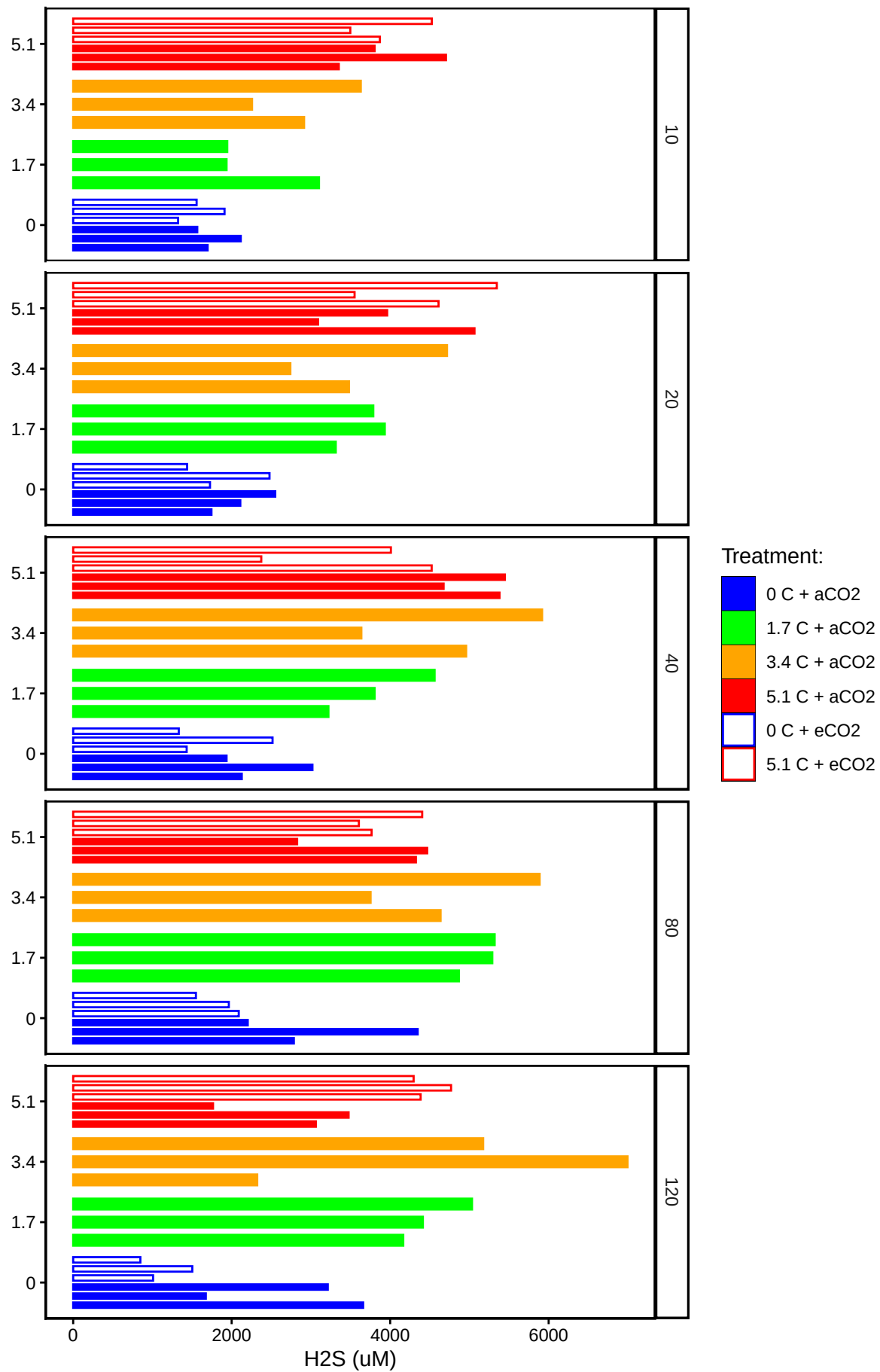


C4: Porewater H2S

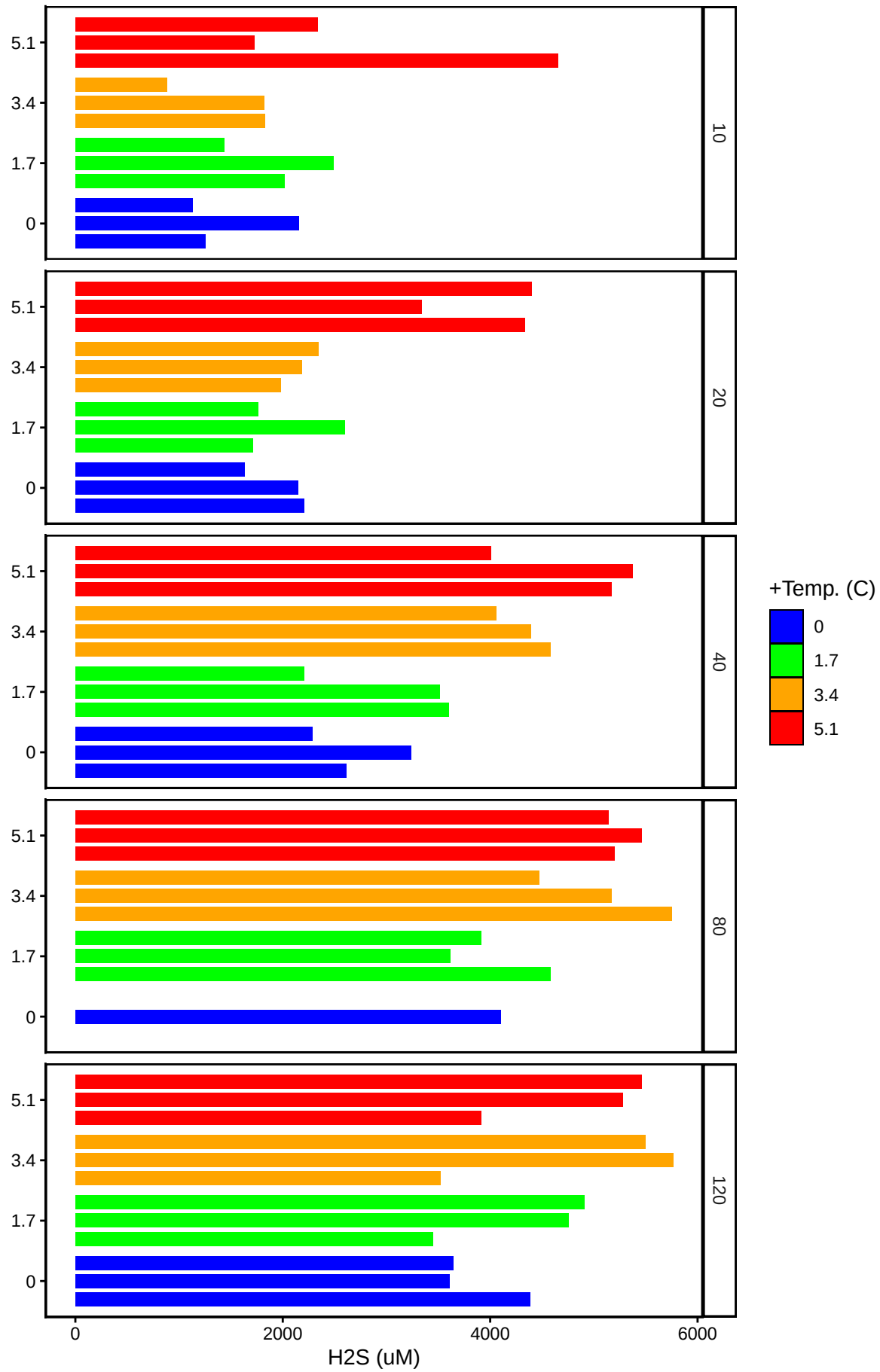


```
## Warning: Using 'size' aesthetic for lines was deprecated in ggplot2 3.4.0.
## i Please use 'linewidth' instead.
## This warning is displayed once per session.
## Call 'lifecycle::last_lifecycle_warnings()' to see where this warning was
## generated.
```

C3: Porewater H₂S by Replicate

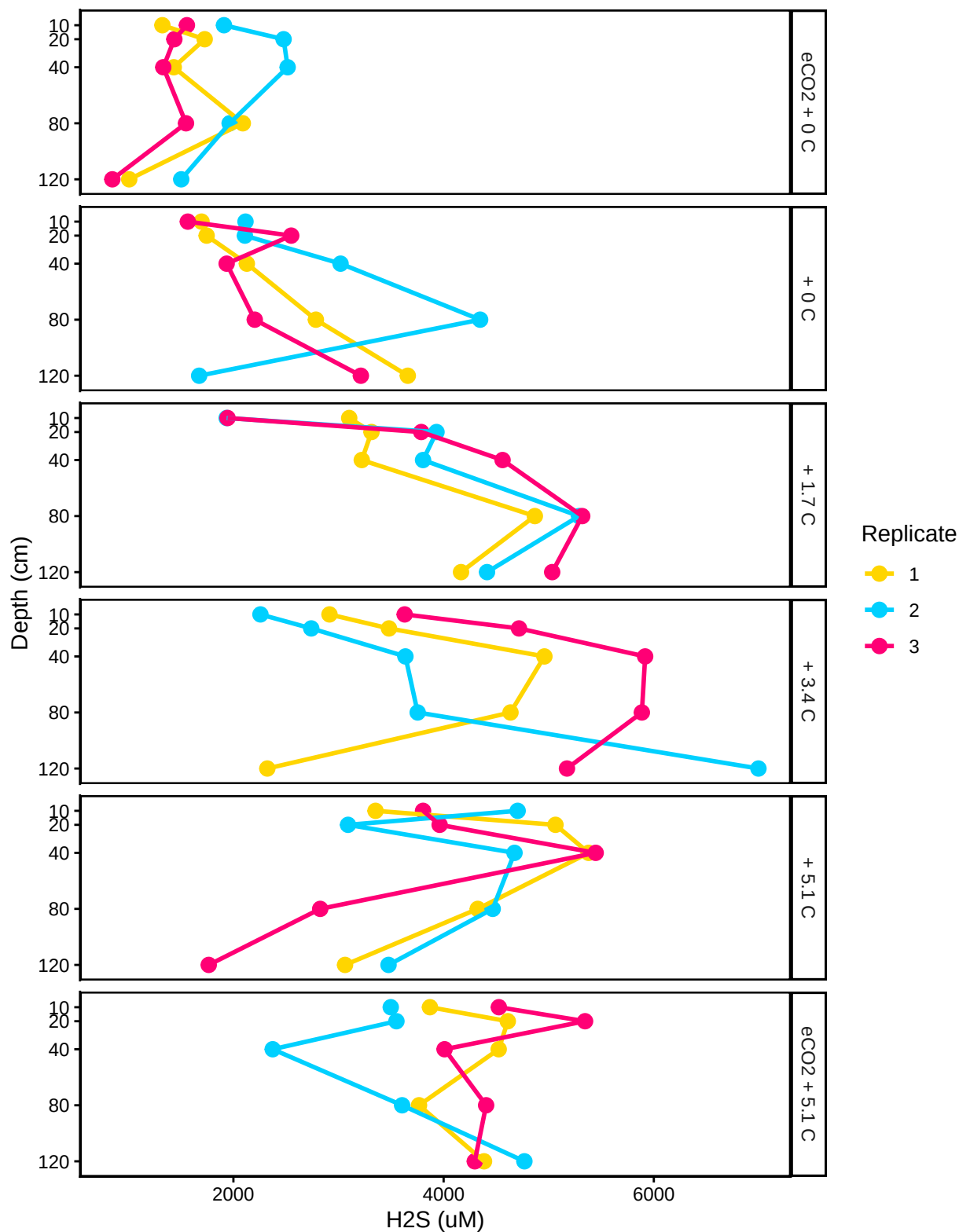

C4: Porewater H2S by Replicate

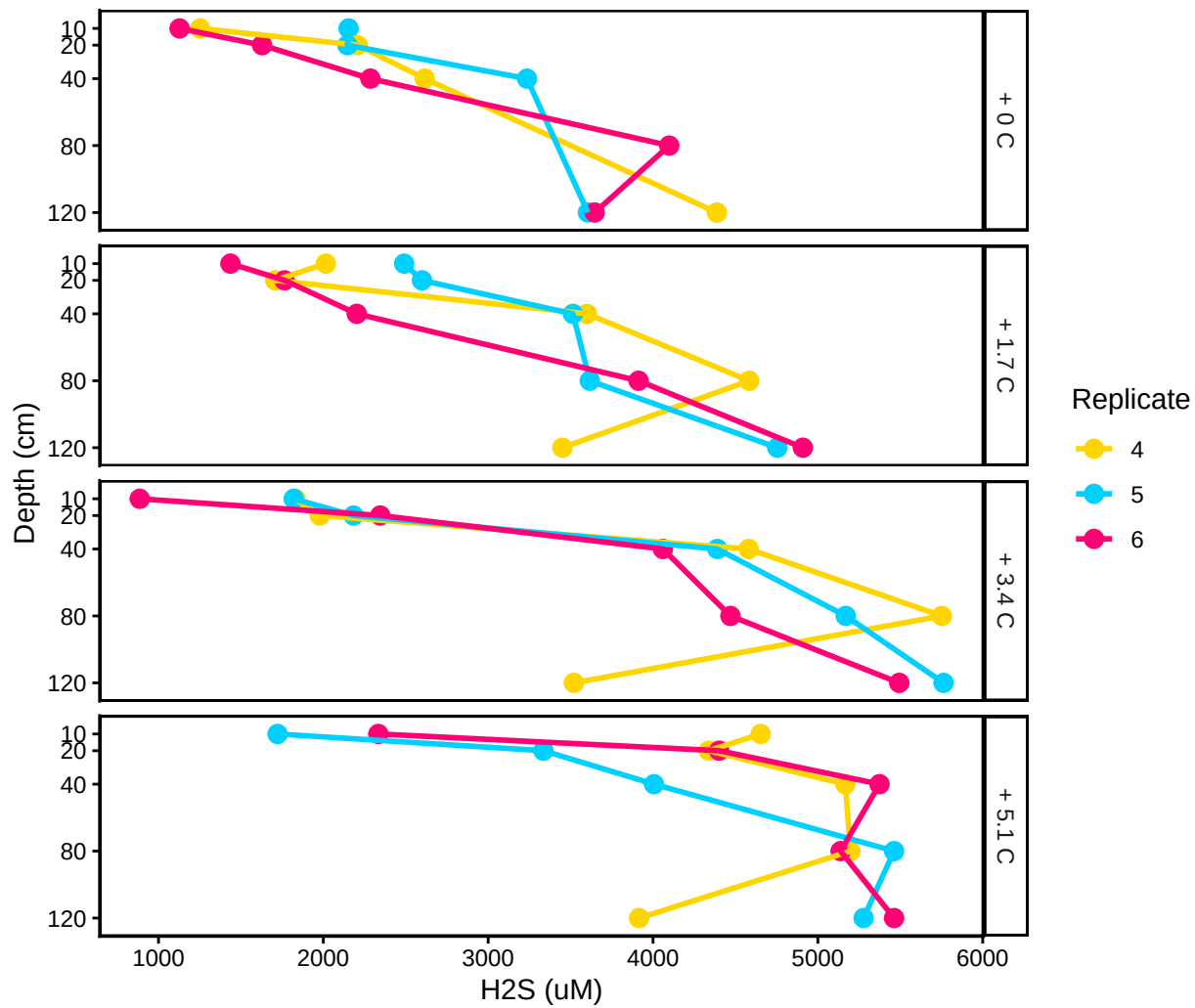

##Export Data

Table 7: Reruns

Sample_ID.x	Conc	Conc_Flag	CV_Flag	Dilution	CV
310-40	2127.645	Within_Range	CV_High	100	11.05918
462-120	5493.686	Within_Range	CV_High	100	13.39841

END