

GENX: Porewater Sulfide

September 2026 Samples

2026-09-29

Code Set up

run information

```
##Before or After Reruns (YES or NO)
After_Reruns = "YES"
#If NO will not write final data file
#And Will Write a Rerun File

###things that need to be changed

sample_month = "September"
sample_year = "2026"
Run_Date = as.Date("9/17/2026", format = "%m/%d/%Y")
Run_by = "Melanie Giessner" #Instrument user
Script_run_by = "Melanie Giessner" #Code user

Run_notes="Ch11-80 not collected" #any notes from run

#Import Plate Data
folder_path_rawdata <- paste0("Raw Data")
Raw_plates <- list.files(path = folder_path_rawdata, full.names = TRUE, pattern = ("(202609|202610)(.+Plate
print(Raw_plates)

## [1] "Raw Data/20260917_LTREB_H2S_Plate1.csv"
## [2] "Raw Data/20260917_LTREB_H2S_Plate2.csv"
## [3] "Raw Data/20260917_LTREB_H2S_Plate3.csv"
## [4] "Raw Data/20260917_LTREB_H2S_Plate4.csv"
## [5] "Raw Data/20260917_LTREB_H2S_Plate5.csv"
## [6] "Raw Data/20260918_LTREB_H2S_Plate6.csv"
## [7] "Raw Data/20260918_LTREB_H2S_Plate7.csv"
## [8] "Raw Data/20260918_LTREB_H2S_Plate8.csv"
## [9] "Raw Data/20260929_LTREB_SMARTX_H2S_Rerun1.csv"

# Define the file path for QAQC log file - NO Need to change just check year
log_path <- here("Porewater QAQC Logs", "GCRew_H2S_QAQC_Log.csv")

# Define final path for data be sure to change year and month and run number
final_path <- "Processed Data/GCRew_GENX_Porewater_H2S_202609.csv"
rerun_path <- "Processed Data/GCRew_GENX_Porewater_H2S_202609_reruns.csv"



##Read in metadata and create similar sample IDs for matching to samples
##Read in data
```


```
##Remove High CV Stds
```

```
##Check R2 & Remove Standards if Necessary
```

```
## Joining with 'by = join_by(Plate, IDs)'
```

Plot Standard Curves

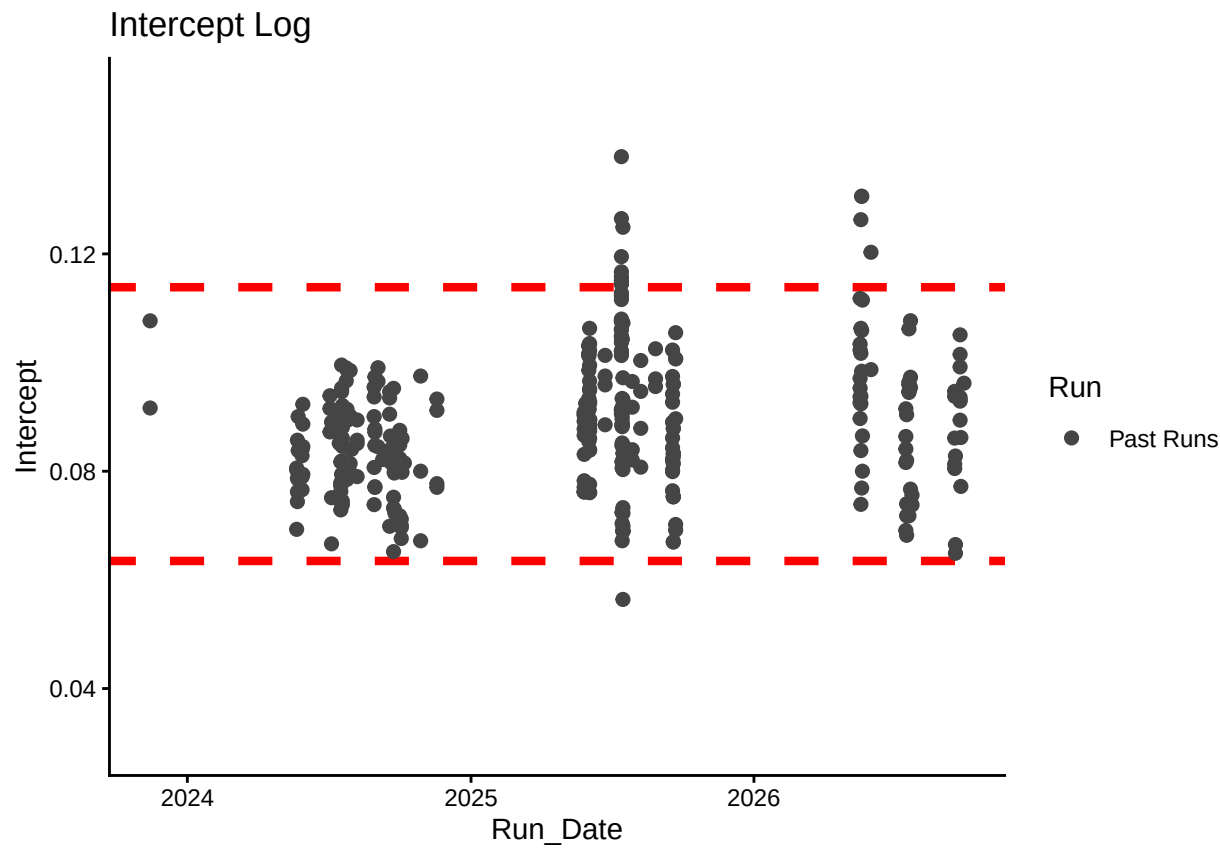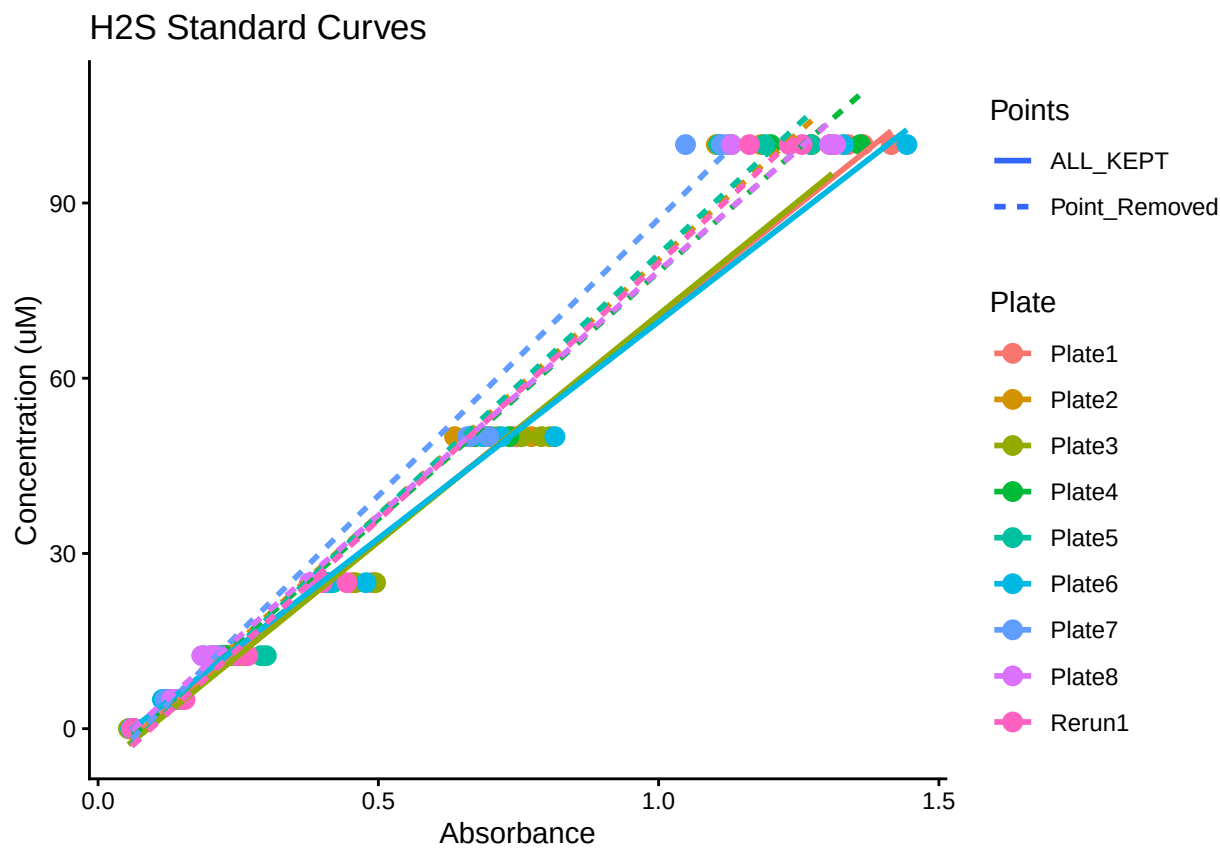

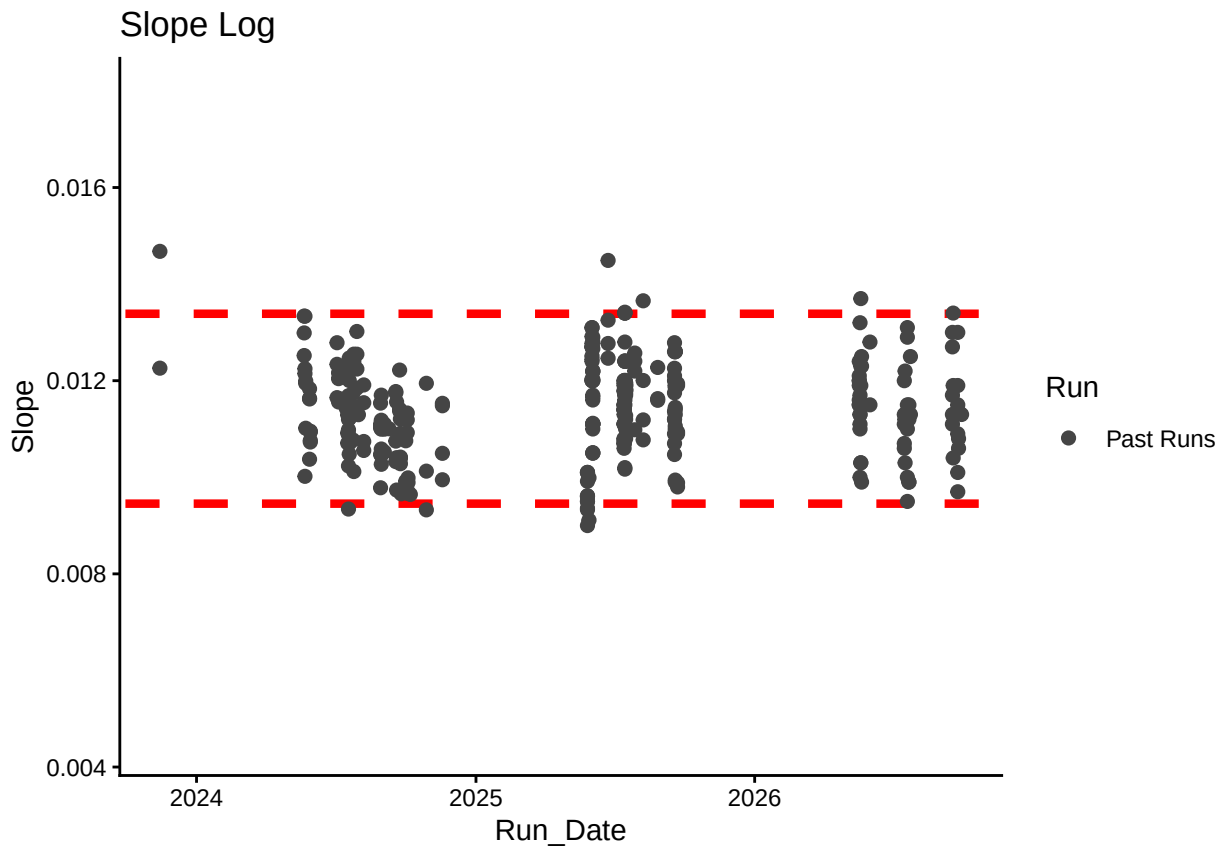


Return Flags For Standard Curve


Plate	R2_Flag	CV_Flag
Plate1	R2_Good	CV_Good
Plate2	R2_Bad	CV_Good
Plate2	R2_Bad	CV_High
Plate3	R2_Good	CV_Good
Plate4	R2_Good	CV_Good
Plate5	R2_Good	CV_Good
Plate6	R2_Good	CV_Good
Plate7	R2_Good	CV_Good
Plate8	R2_Good	CV_Good
Rerun1	R2_Good	CV_Good

Plate	Curve_RMVD
Plate1	KEPT
Plate2	REMOVED
Plate3	KEPT
Plate4	KEPT
Plate5	KEPT
Plate6	REMOVED
Plate7	KEPT
Plate8	KEPT
Rerun1	KEPT

Method Detection Limit

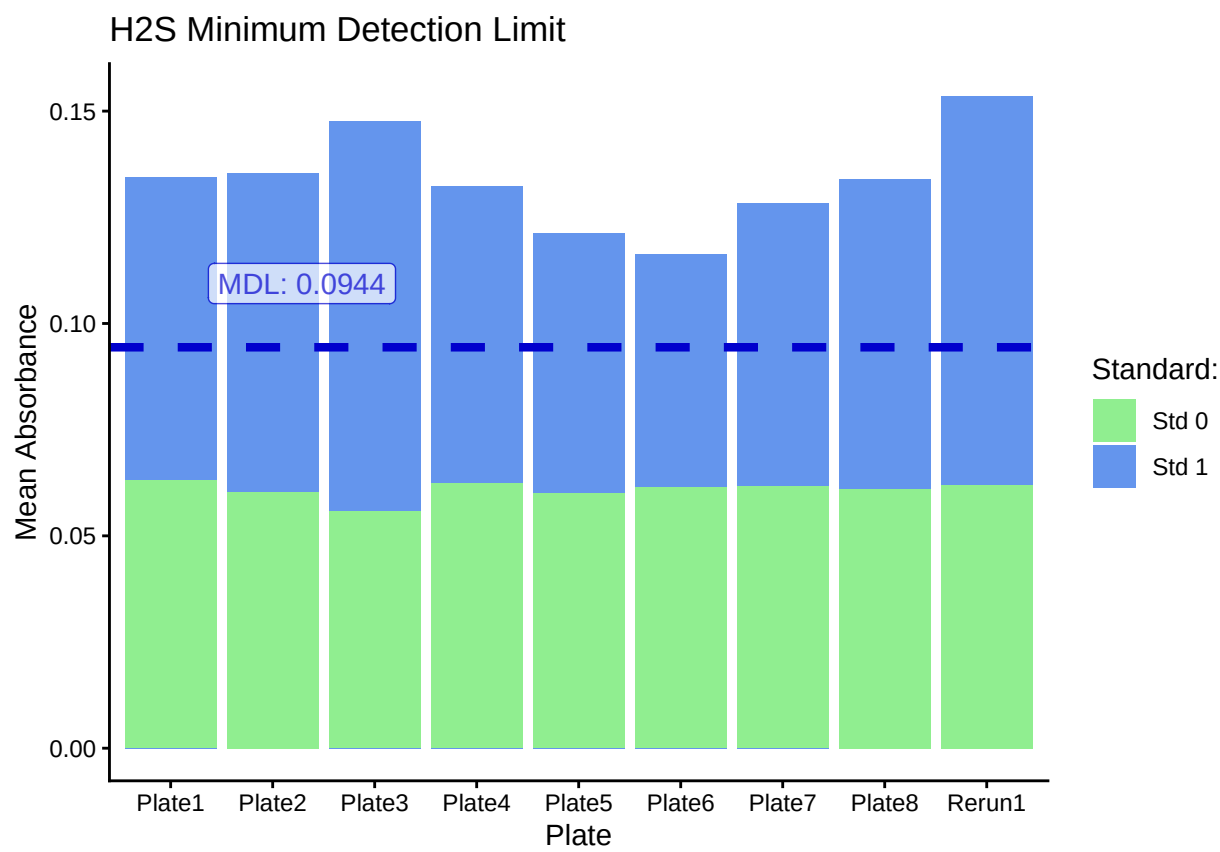

Compare Check Standards to Standards

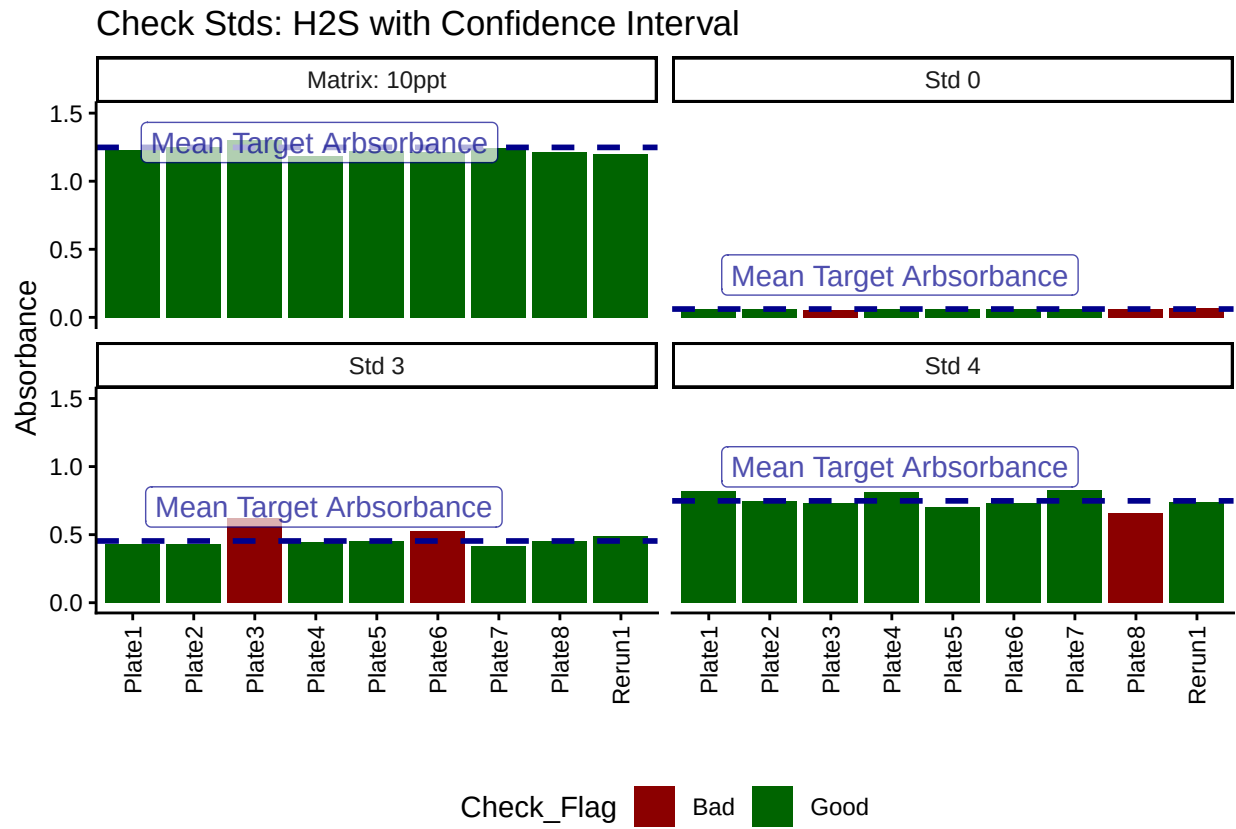



Display Any Check Flags

Plate	CHKs_Flag
Plate3	Std0Chk_Bad
Plate3	Std3Chk_Bad
Plate6	Std3Chk_Bad
Plate8	Std0Chk_Bad
Plate8	Std4Chk_Bad
Rerun1	Std0Chk_Bad


Plate	CHK_CV_Flags
Plate8	Std3Chk_CV_High
Rerun1	Std0Chk_CV_High

```
##Calculate Sulfide Concentrations & Add Concentration Flags
##Remove High CV Samples and Average
```

Plot Samples

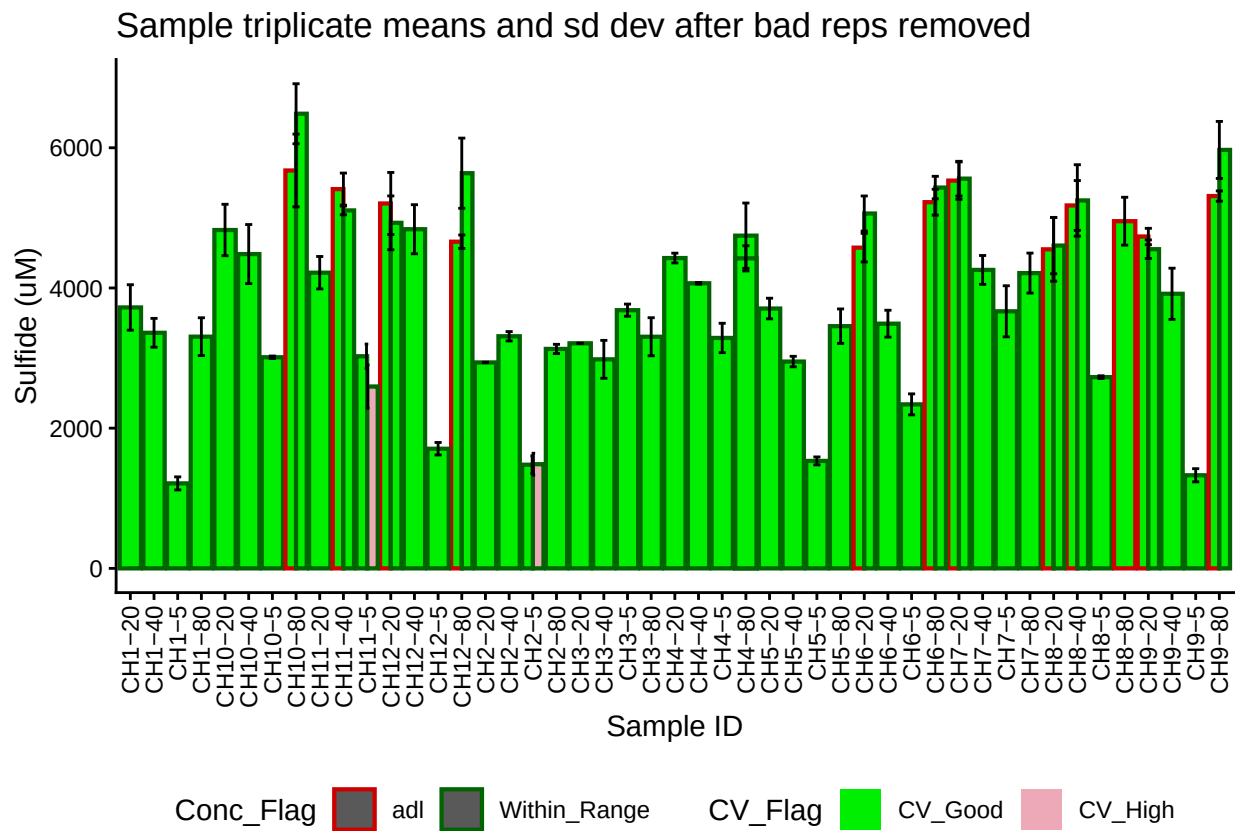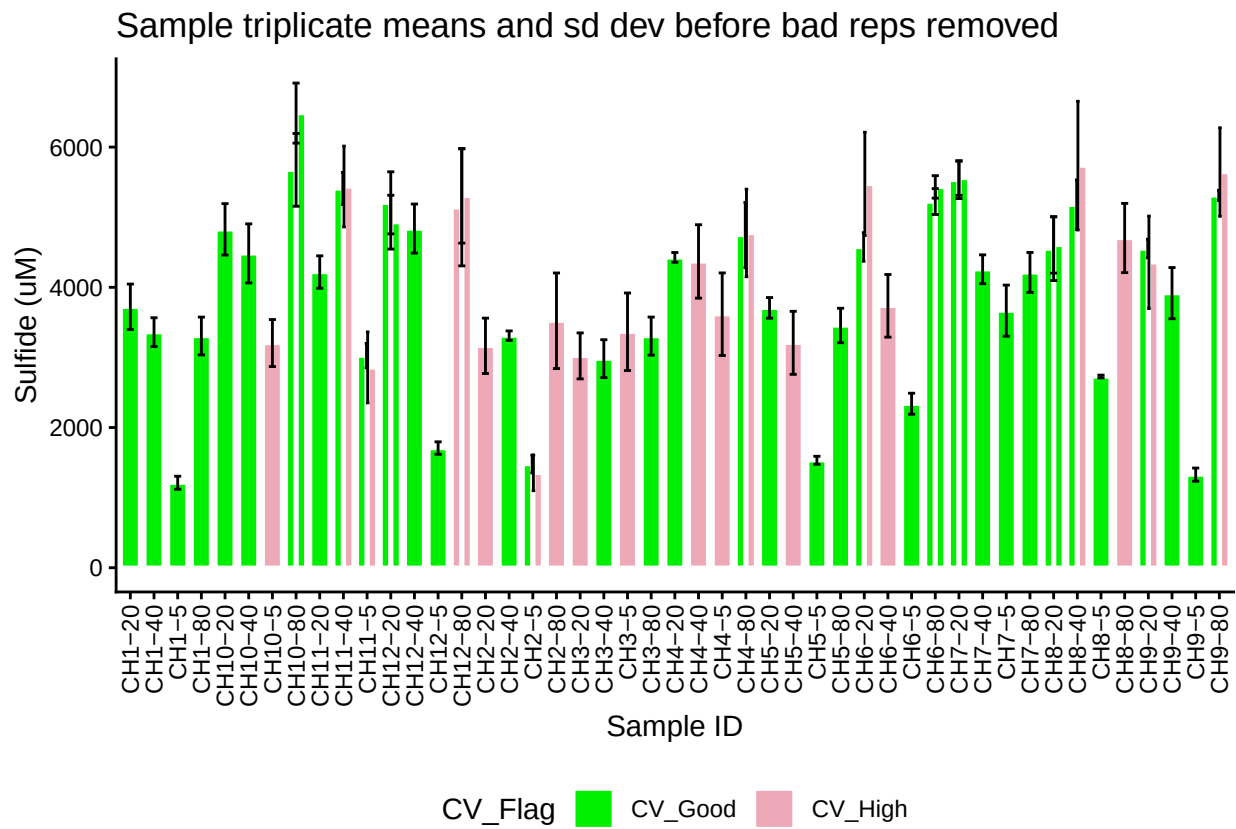

Duplicate QAQC

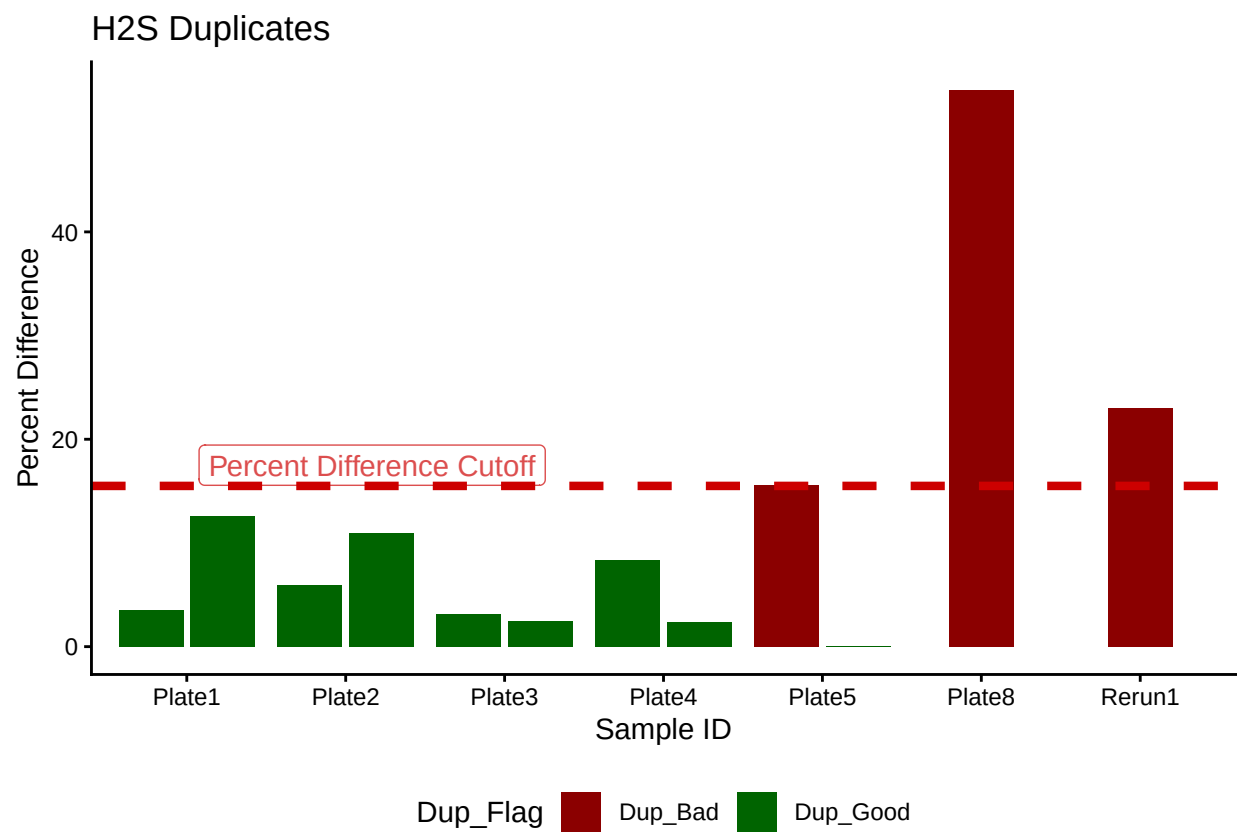


Display Dup Flags

[1] ">60% Dups Pass"

Plate	Dup_Flag
Plate1	Dup_Good
Plate2	Dup_Good
Plate3	Dup_Good
Plate4	Dup_Good
Plate5	Dup_Bad
Plate8	Dup_Bad
Rerun1	Dup_Bad

Spike QAQC

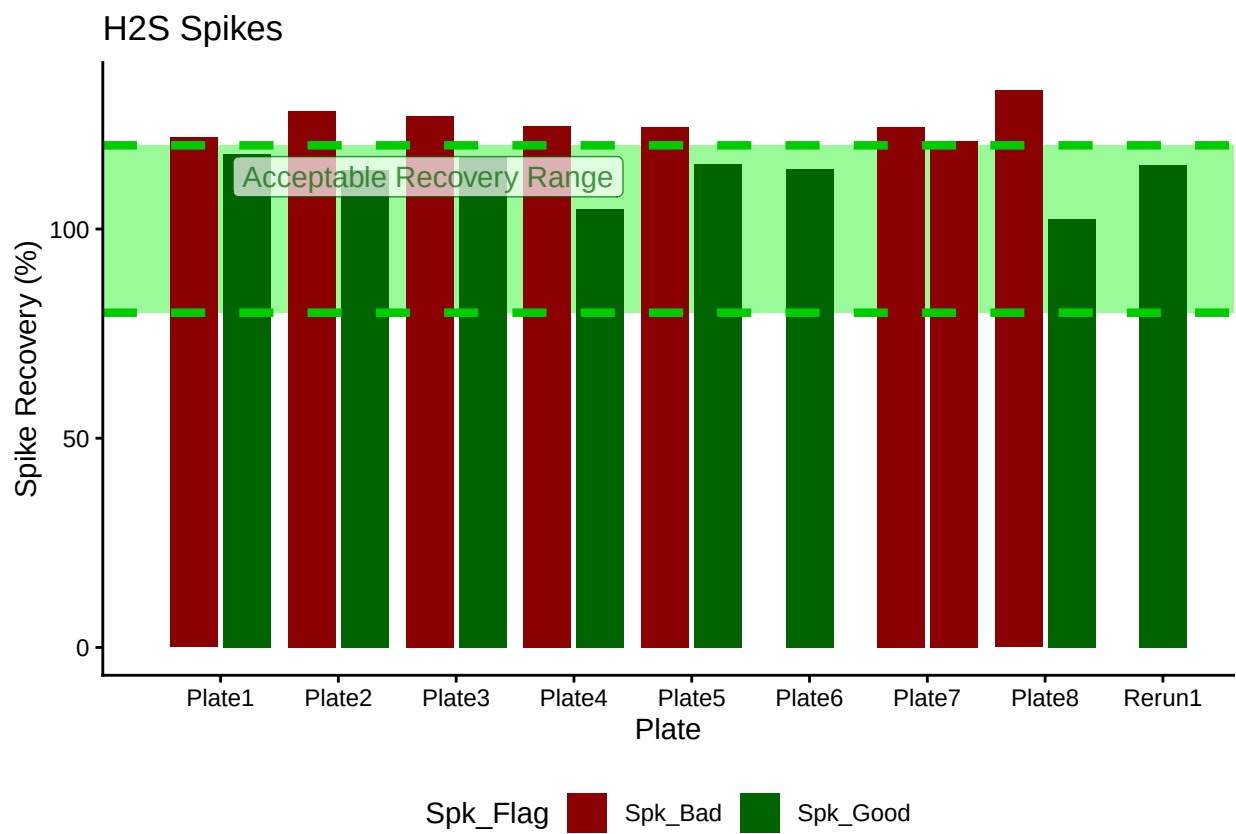

Display Spike Flags

```
## [1] "Spike Missing Pair"
```

Plate	Spike_Flag
Plate1	Spk_Bad
Plate2	Spk_Bad
Plate3	Spk_Bad
Plate4	Spk_Bad
Plate5	Spk_Bad
Plate6	Spk_Good
Plate7	Spk_Bad
Plate8	Spk_Bad
Rerun1	Spk_Good

```
## [1] "<60% of Spikes Pass, Rerun"
```

Merge samples with Metadata & check all samples are present

```
## Some sample IDs are missing from data.
```

```
## [1] "Ch11-80"
```

```
##Add Flags
```


```
##Organize Data
```

```
## Adding missing grouping variables: 'Sample_ID.x'  
## Adding missing grouping variables: 'Sample_ID.x'
```


Visualize Data

```
## Visualize Data
```

GENX: Porewater H2S

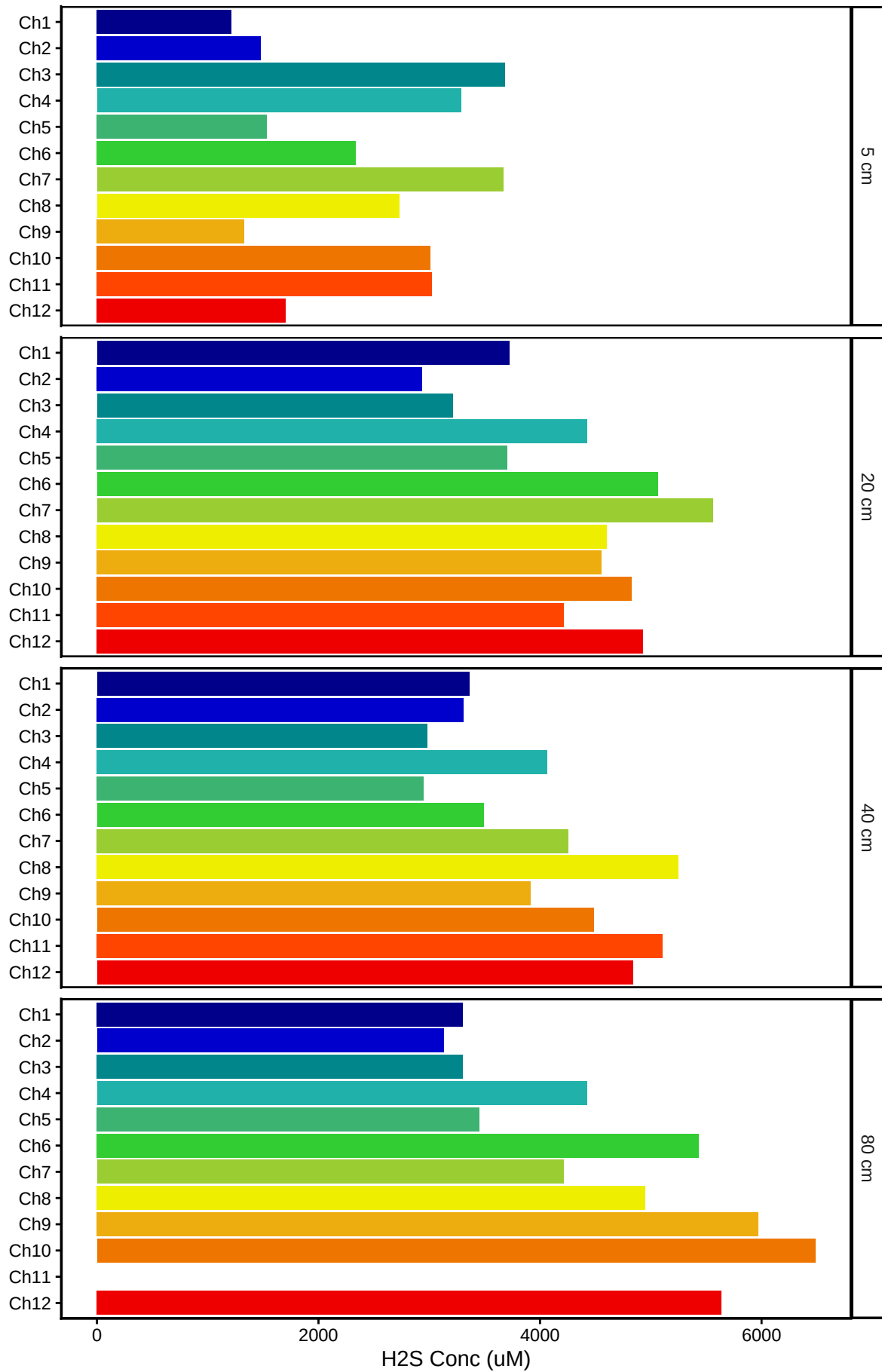


GENX: Porewater H2S by Chamber

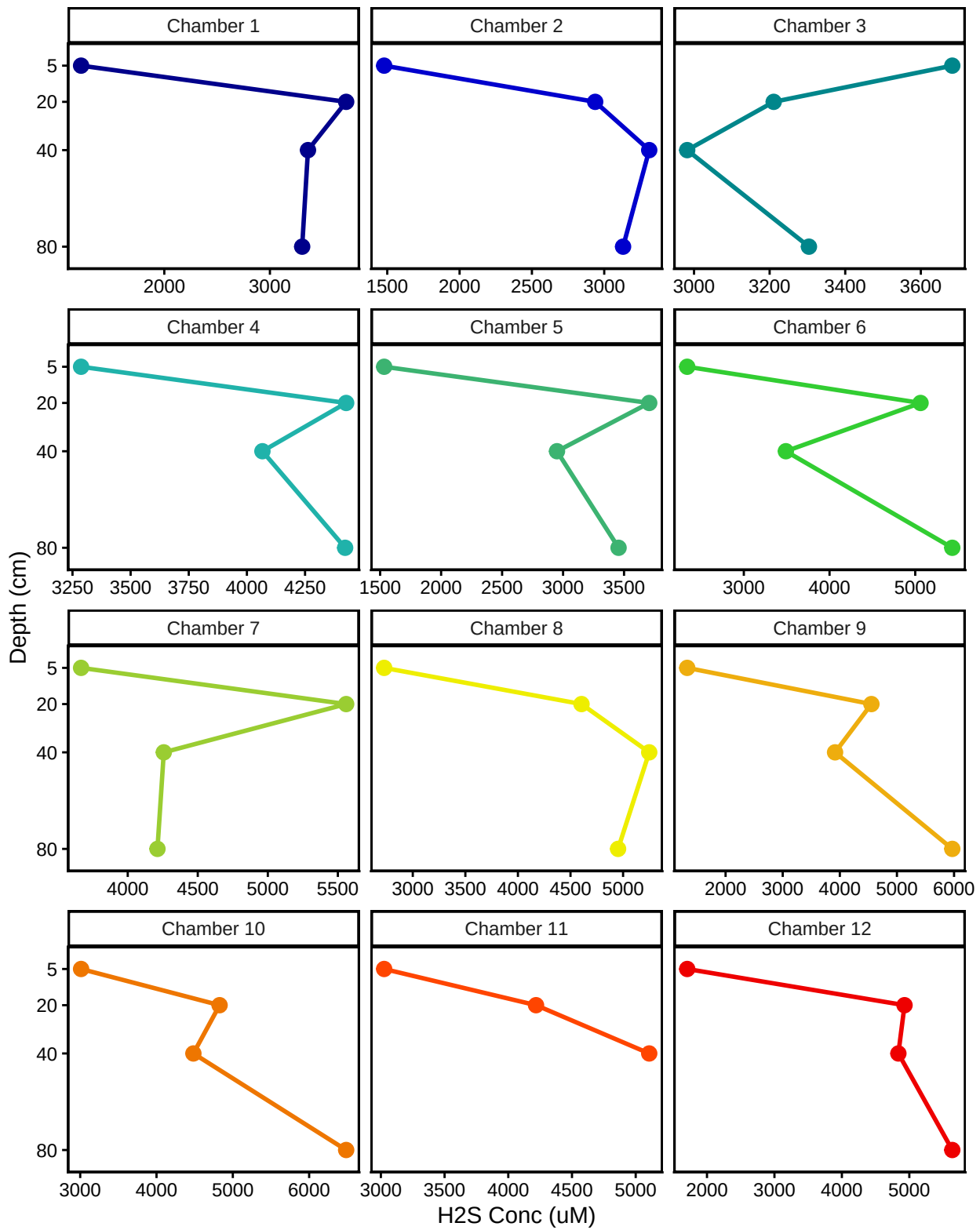


Visualize Data

Visualize Data

GENX: Porewater H2S

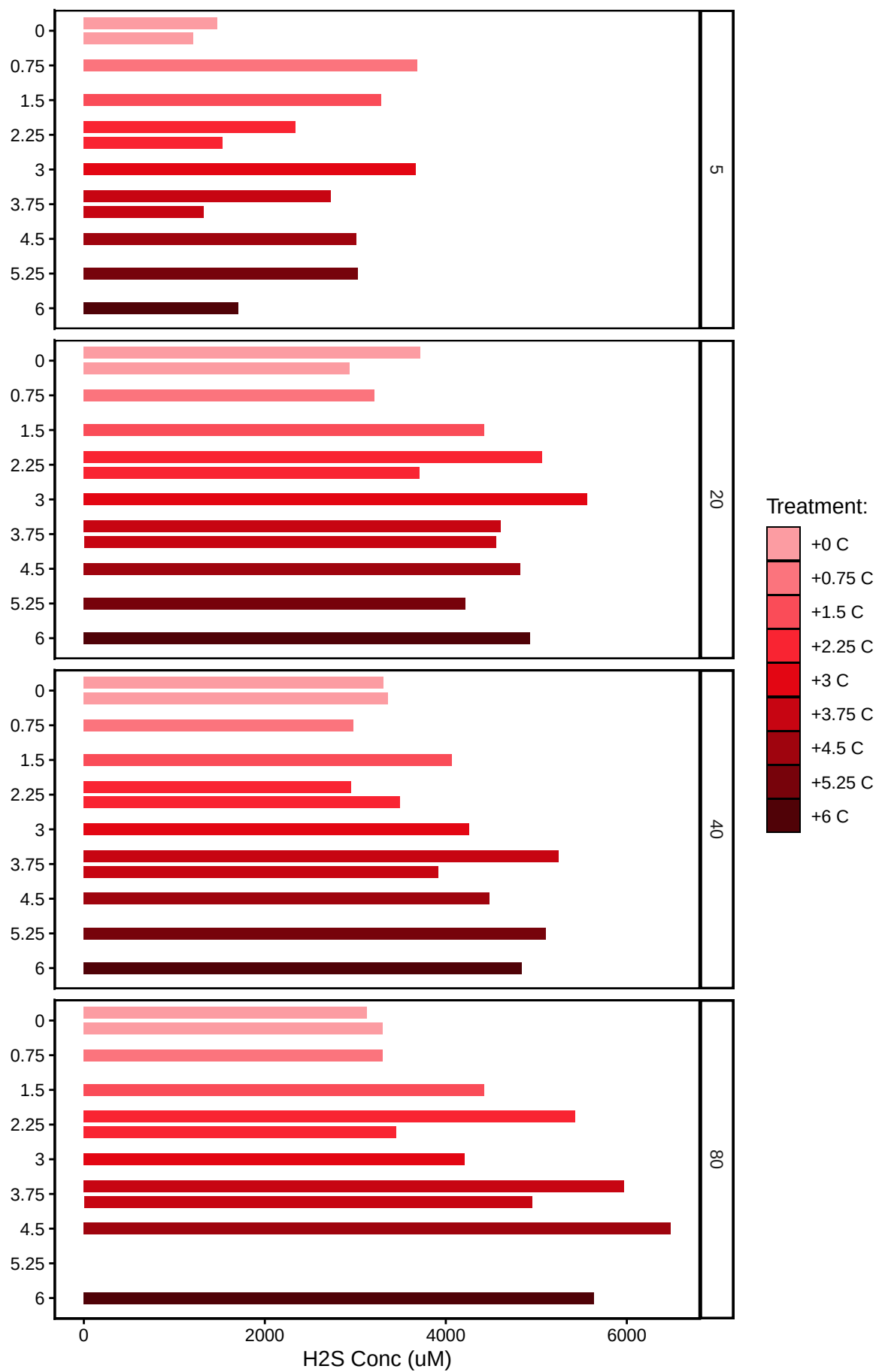

GENX: Porewater H2S by Replicate

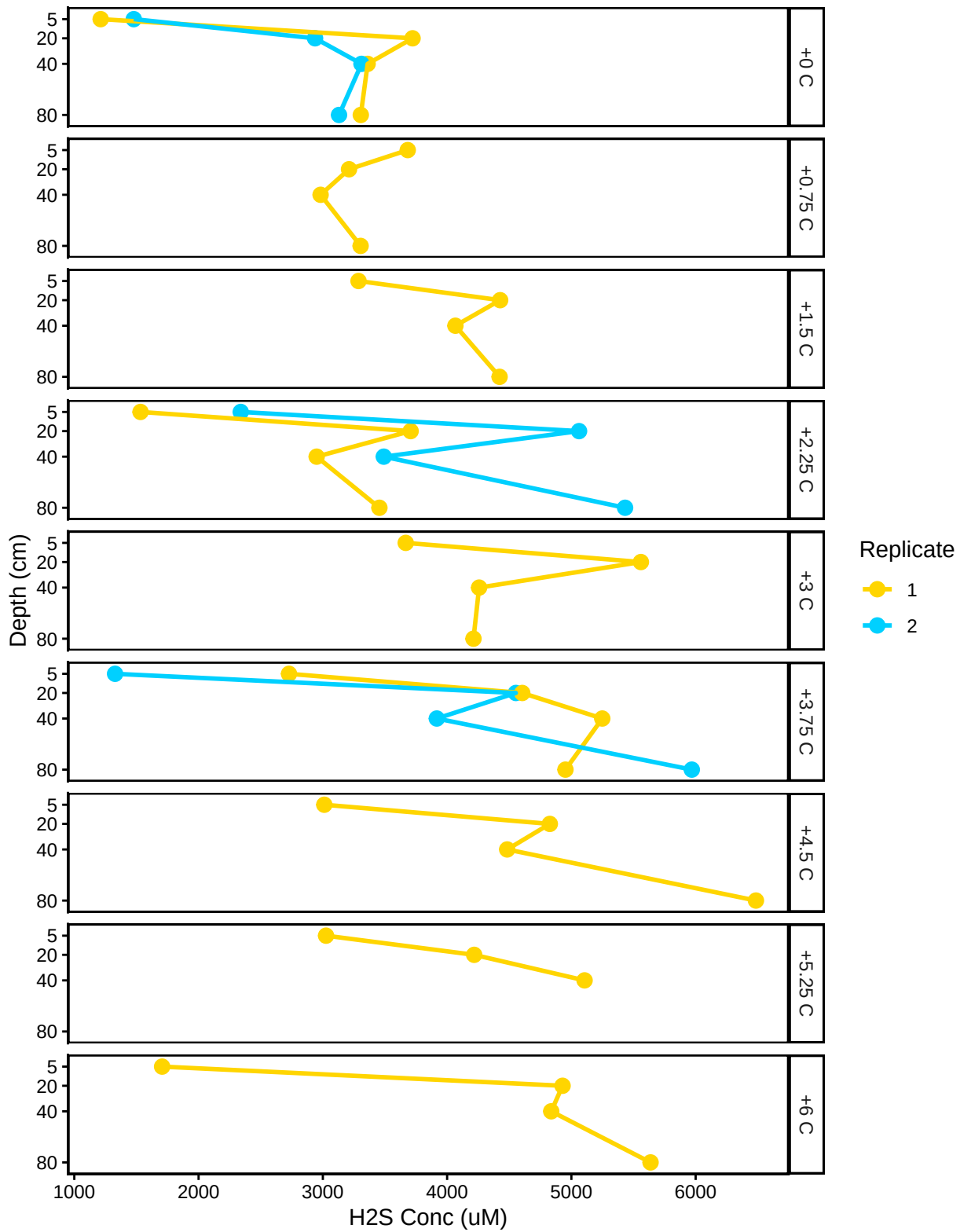

##Export Data

END