

GENX: Porewater Sulfide

Month Year Samples

2026-05-28

Code Set up

run information

```
##Before or After Reruns (YES or NO)
After_Reruns = "NO"
#If NO will not write final data file
#And Will Write a Rerun File

###things that need to be changed

sample_month = "May"
sample_year = "2026"
Run_Date = as.Date("5/20/2026", format = "%m/%d/%Y")
Run_by = "Melanie Giessner" #Instrument user
Script_run_by = "Melanie Giessner" #Code user

Run_notes="Ch1-80 not collected" #any notes from run

#Import Plate Data

plates<- c("Plate1","Plate2","Plate3","Plate4")

folder_path_rawdata <- paste0("Raw Data")
Raw_plates <- list.files(path = folder_path_rawdata, full.names = TRUE,pattern = "Plate")
head(Raw_plates)

## [1] "Raw Data/20260519_GENX_H2S_Plate1.csv"
## [2] "Raw Data/20260519_GENX_H2S_Plate2.csv"
## [3] "Raw Data/20260519_GENX_H2S_Plate3.csv"
## [4] "Raw Data/20260520_SMARTX_H2S_Plate6.csv"

#Import Run Date for Each Plate for Std curve qaqc
Run_dates <- tibble(
  Plate = plates,
  Run_Date = c("05/19/2026","05/19/2026","05/19/2026","05/20/2026")
)

Run_dates <- Run_dates %>% mutate(Run_Date = as.Date(Run_Date, format = "%m/%d/%Y"))

# Define the file path for QAQC log file - NO Need to change just check year
log_path <- "Processed Data/SMARTX_H2S_QAQC_2026.csv"
```

```
# Define final path for data be sure to change year and month and run number
final_path <- "Processed Data/GCReW_GENX_Porewater_H2S_202605.csv"
rerun_path <- "Processed Data/GCReW_GENX_Porewater_H2S_202605_reruns.csv"
```

```
##Read in metadata and create similar sample IDs for matching to samples
```



```
##Read in data
```


```
##Remove High CV Stds
```

```
##Check R2 & Remove Standards if Necessary
```

```
## Joining with 'by = join_by(Plate, IDs)'
```

Plot Standard Curves

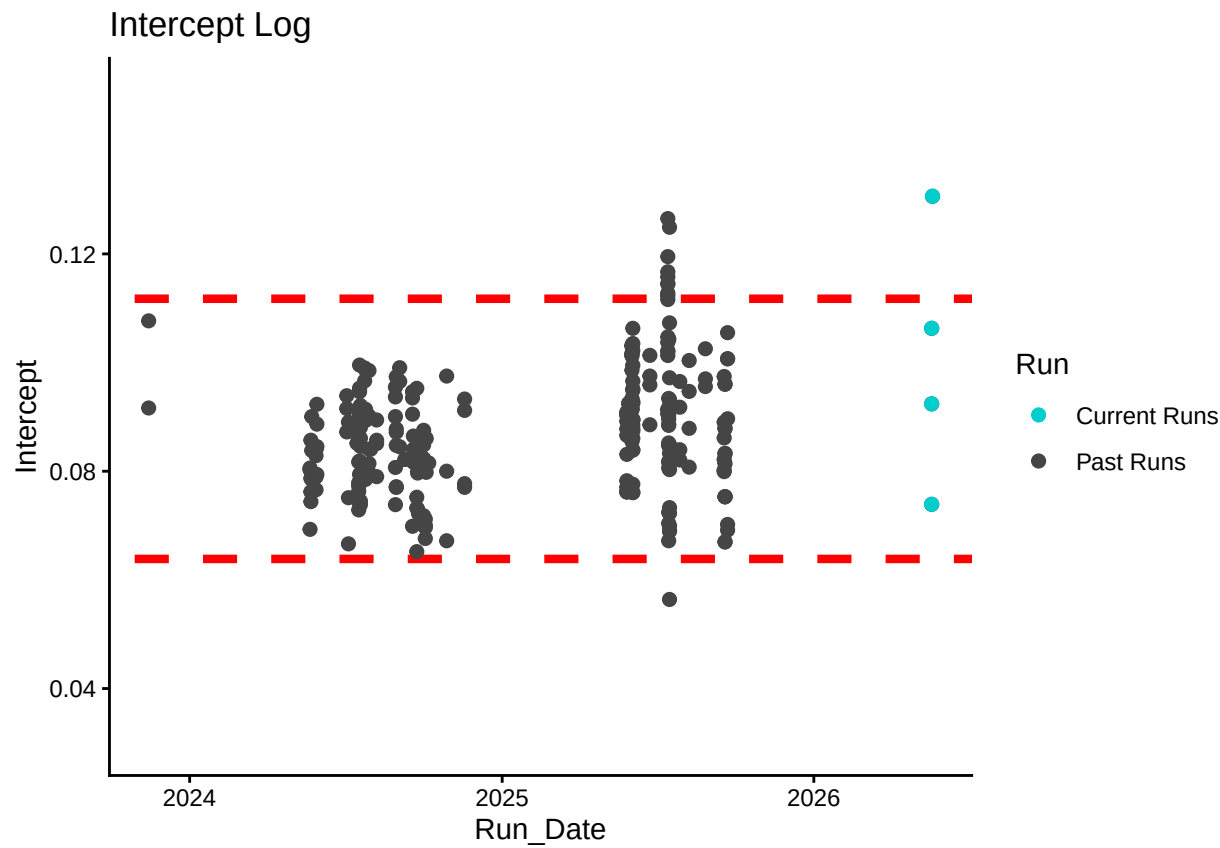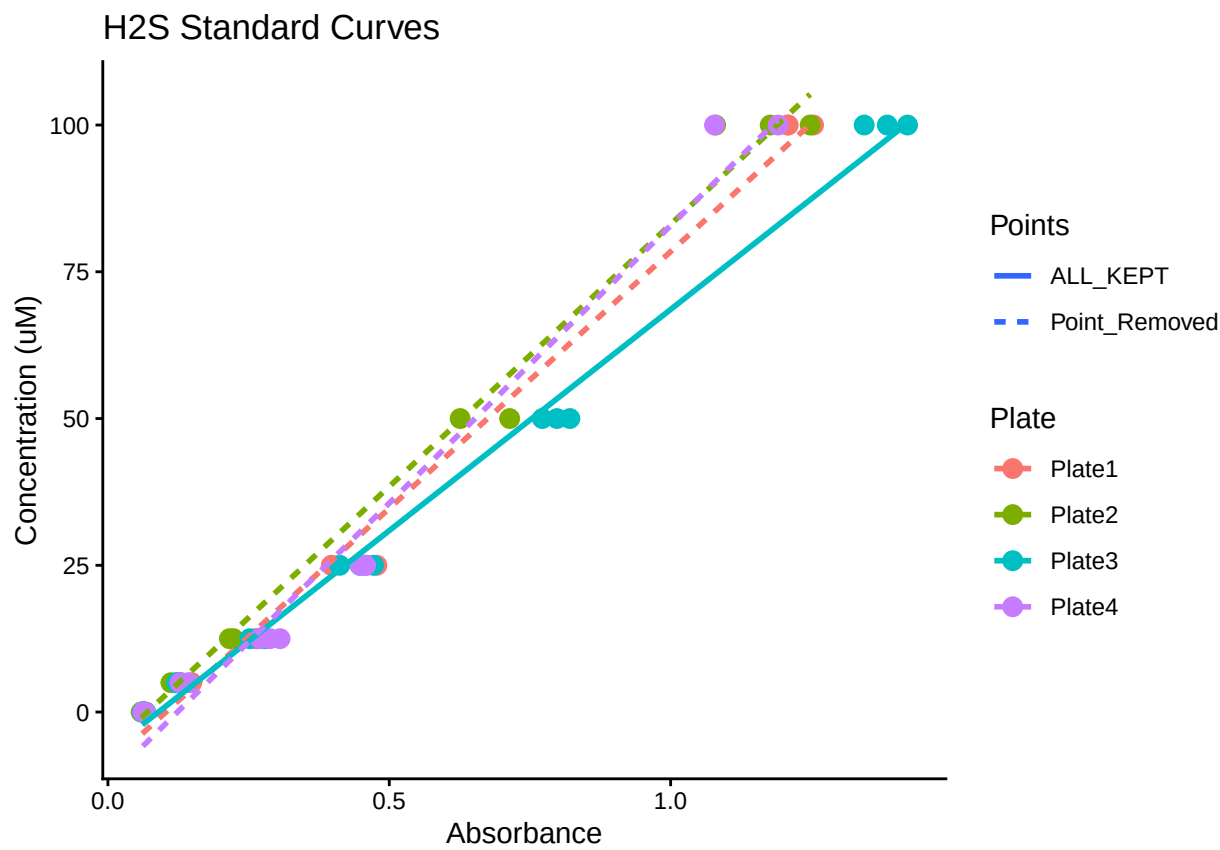

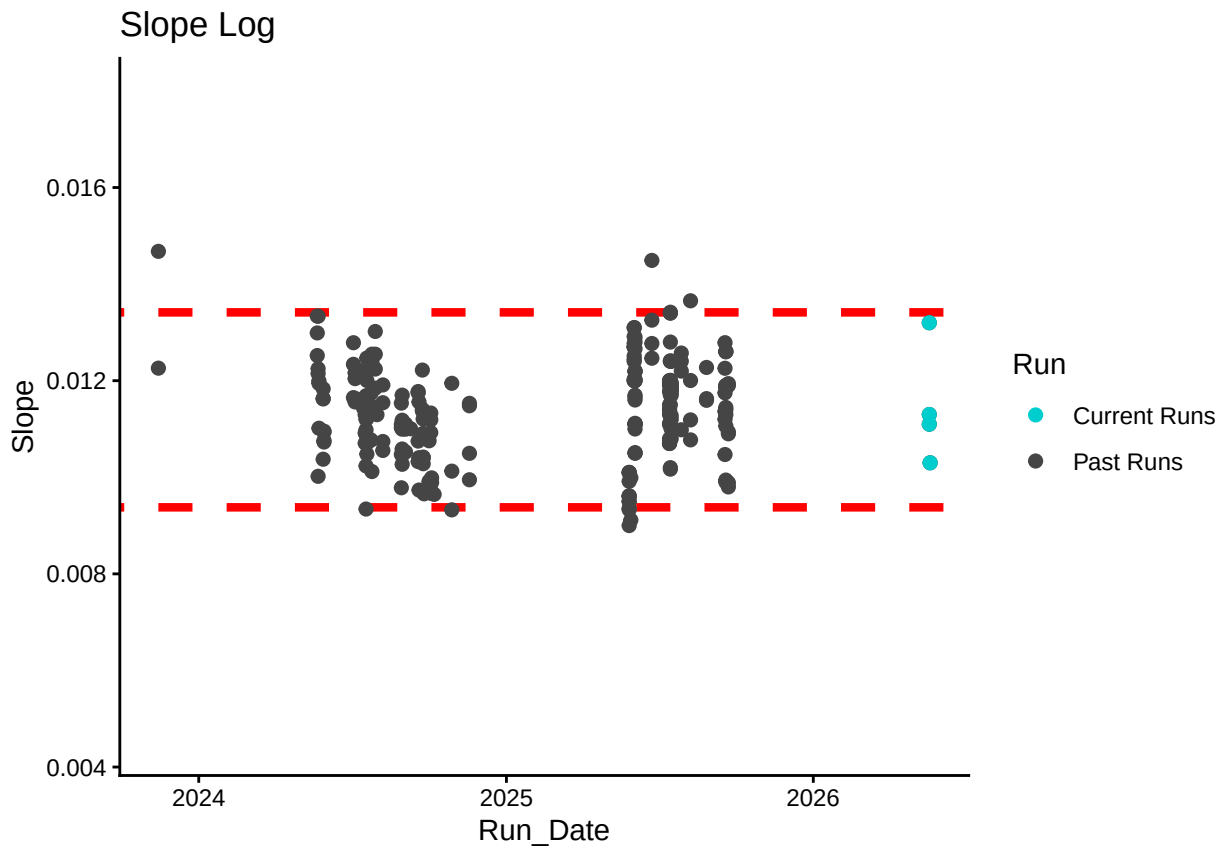


```
## Warning: There was 1 warning in 'mutate()'.
## i In argument: 'Date_Diff = abs(difftime(Date_1, Run_Date, units = "days"))':
## Caused by warning in 'unclass(time1) - unclass(time2)':
## ! longer object length is not a multiple of shorter object length
```

Return Flags For Standard Curve


Plate	R2_Flag	CV_Flag
Plate1	R2_Good	CV_Good
Plate1	R2_Good	CV_High
Plate2	R2_Good	CV_Good
Plate3	R2_Good	CV_Good
Plate4	R2_Bad	CV_Good

Plate	Curve_RMVD
Plate1	KEPT
Plate1	KEPT
Plate2	KEPT
Plate3	KEPT
Plate4	REMOVED

Method Detection Limit

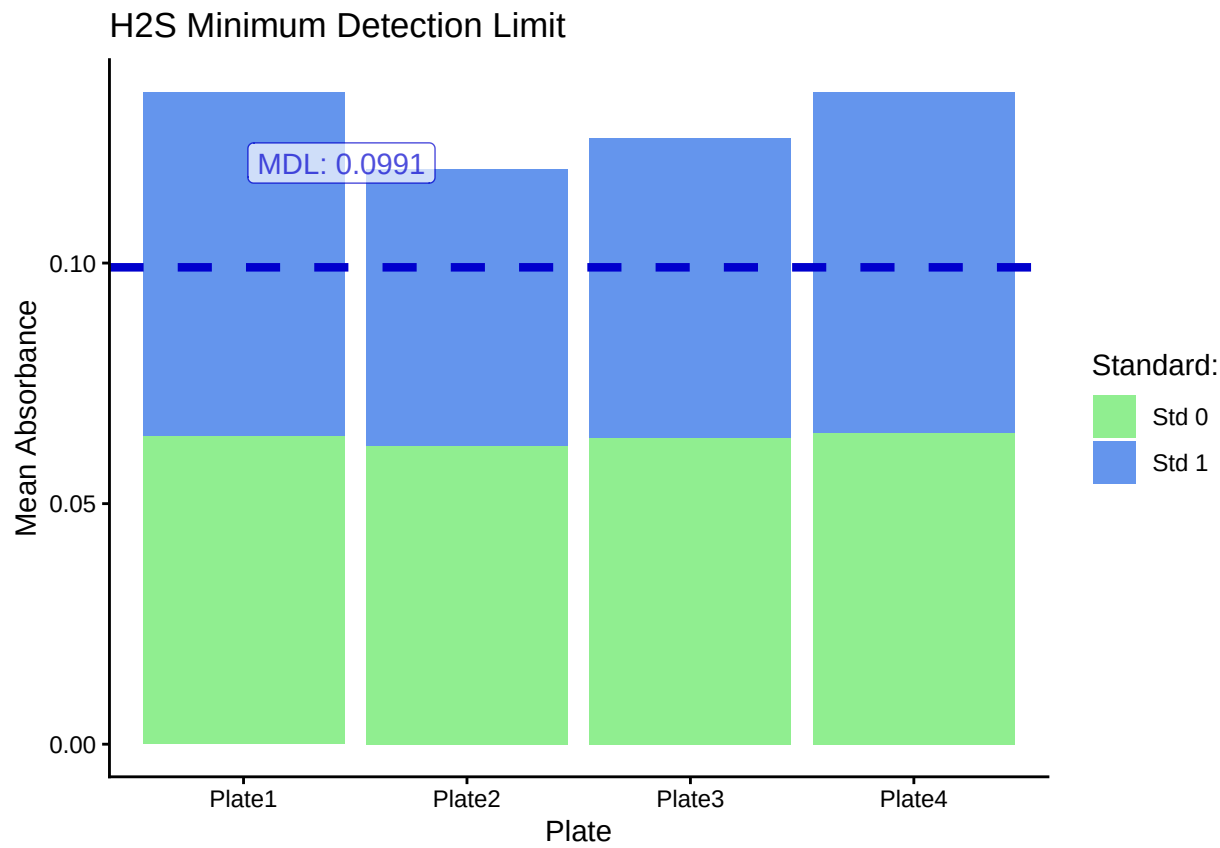

Compare Check Standards to Standards

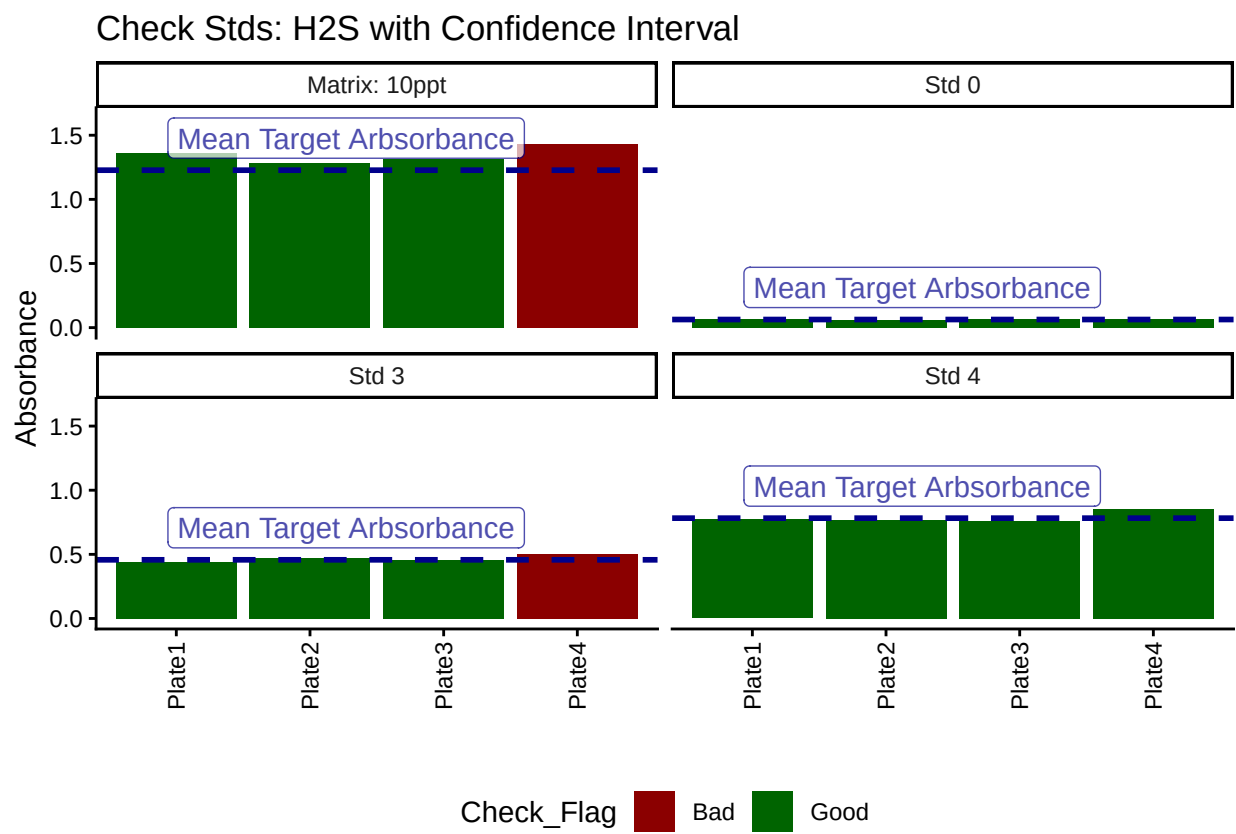

Display Any Check Flags



Plate	CHKs_Flag
Plate4	Std3Chk_Bad
Plate4	Std5Chk_Bad

Plate	CHK_CV_Flags
-------	--------------


##Calculate Sulfide Concentrations & Add Concentration Flags

##Remove High CV Samples and Average

Plot Samples

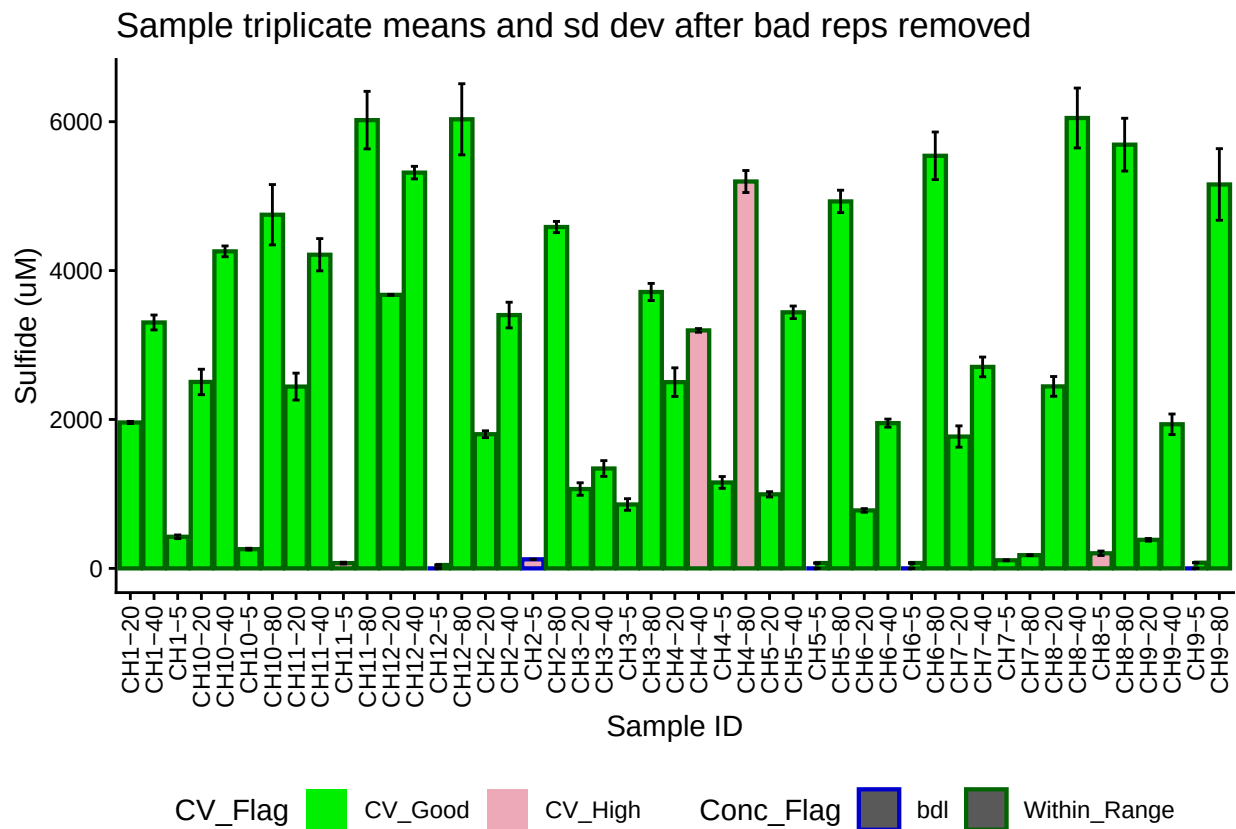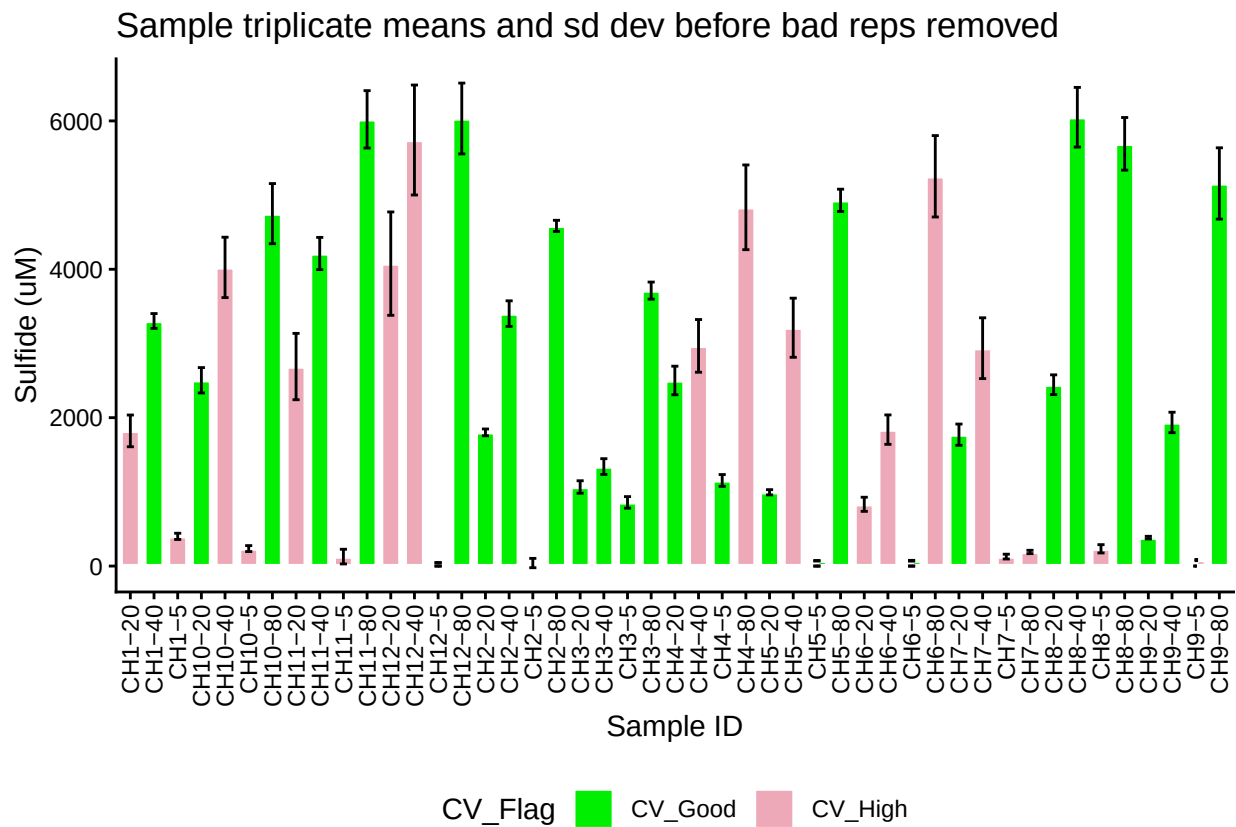

Duplicate QAQC

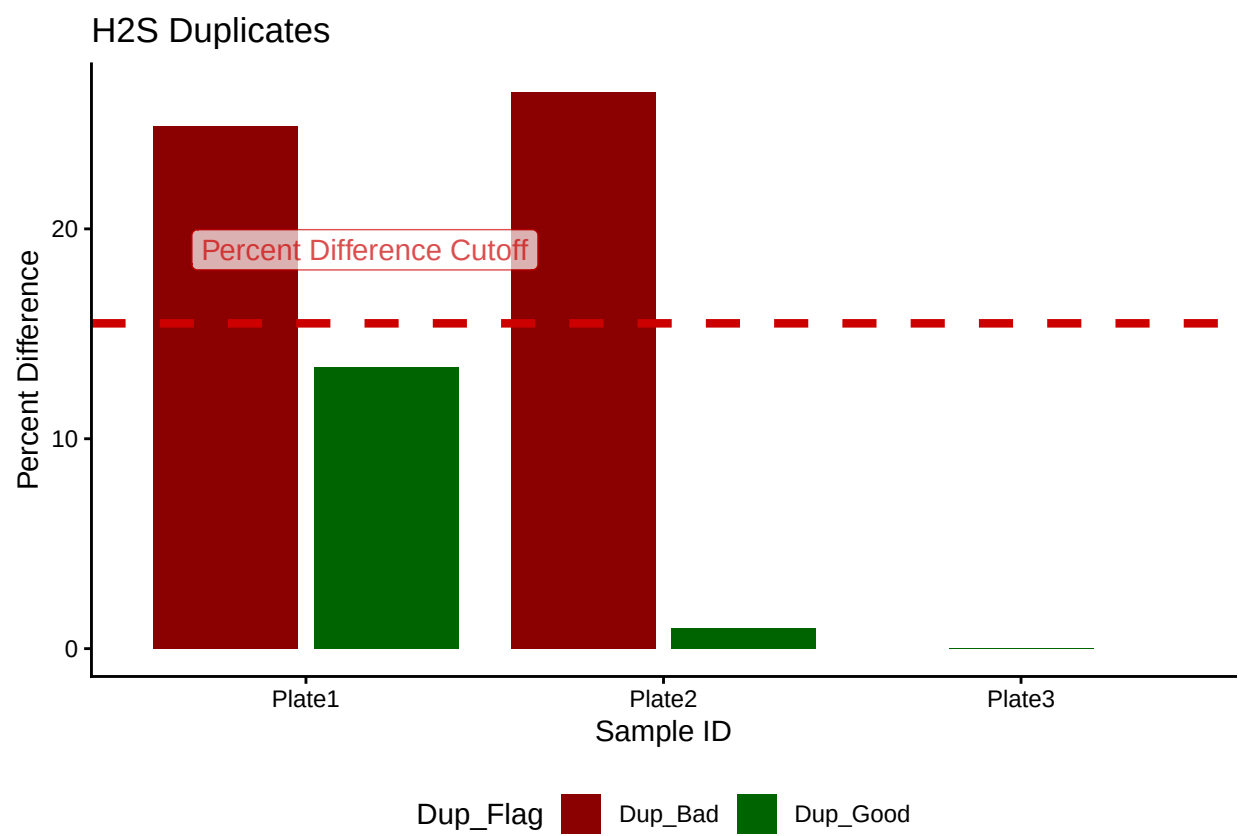


Display Dup Flags

[1] ">60% Dups Pass"

Plate	Dup_Flag
Plate1	Dup_Bad
Plate2	Dup_Bad
Plate3	Dup_Good

Spike QAQC

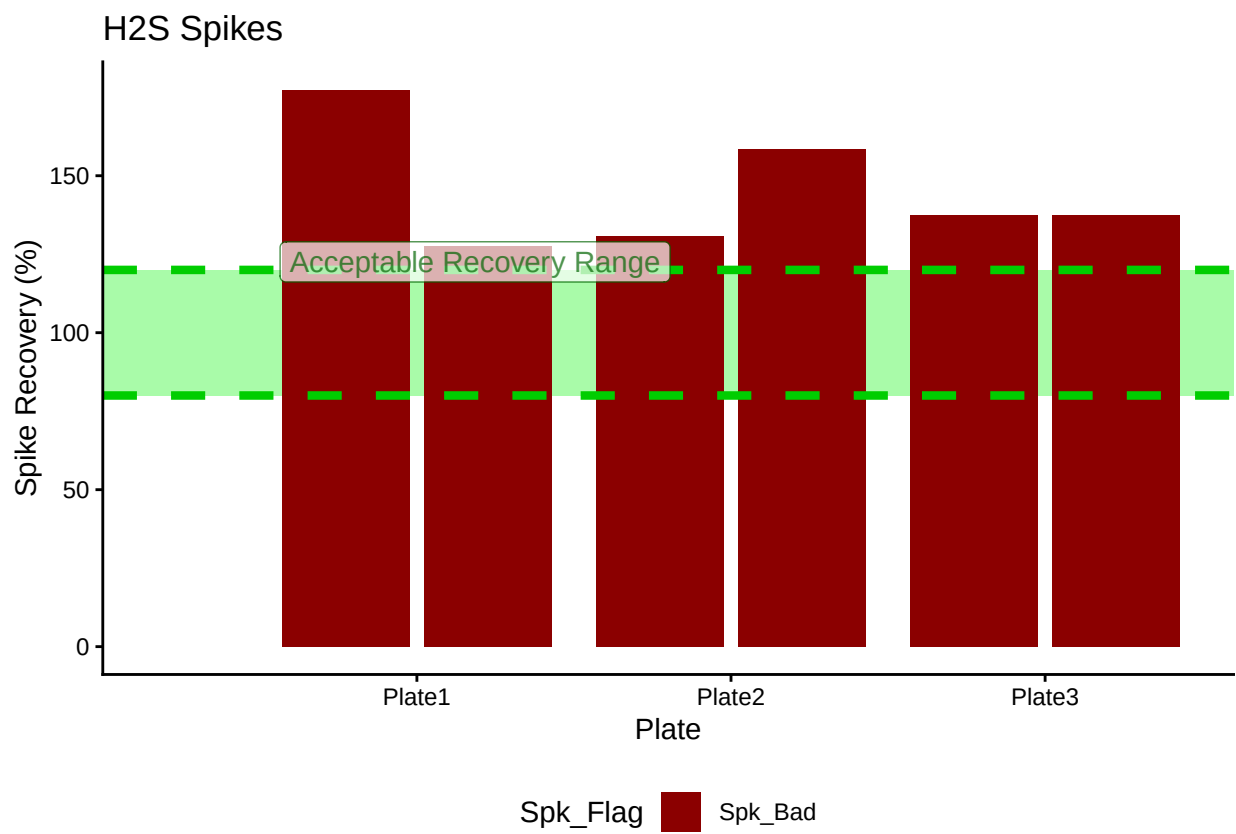

Display Spike Flags

```
## [1] "Spikes Paired"
```

Plate	Spike_Flag
Plate1	Spk_Bad
Plate2	Spk_Bad
Plate3	Spk_Bad

```
## [1] "<60% of Spikes Pass, Rerun"
```

Merge samples with Metadata & check all samples are present

```
## Some sample IDs are missing from data.
```

```
## [1] "Ch1-80"
```


```
##Add Flags
```

```
##Organize Data
```


Visualize Data

```
## Visualize Data
```

GENX: Porewater H2S

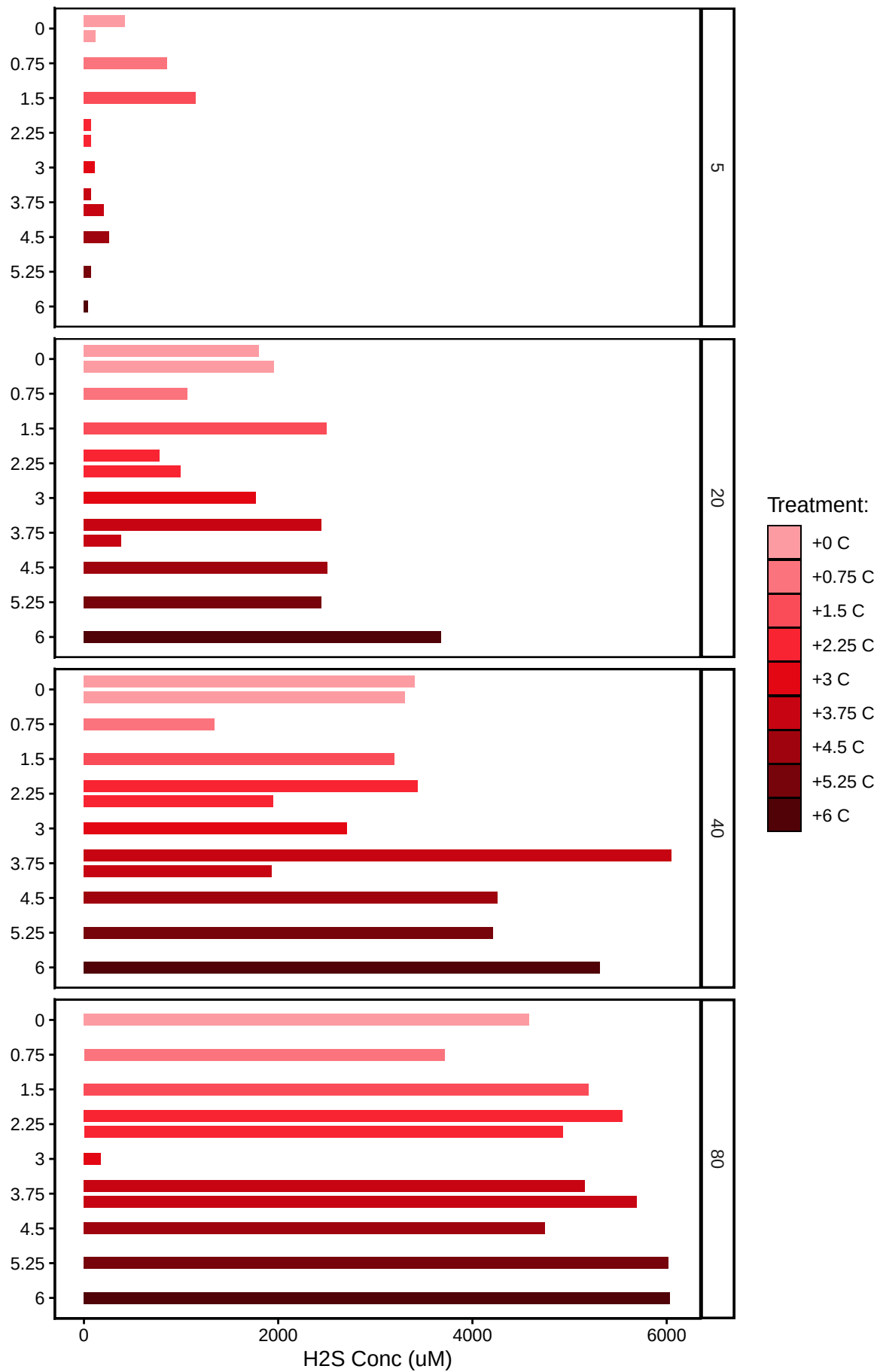

GENX: Porewater H2S by Replicate

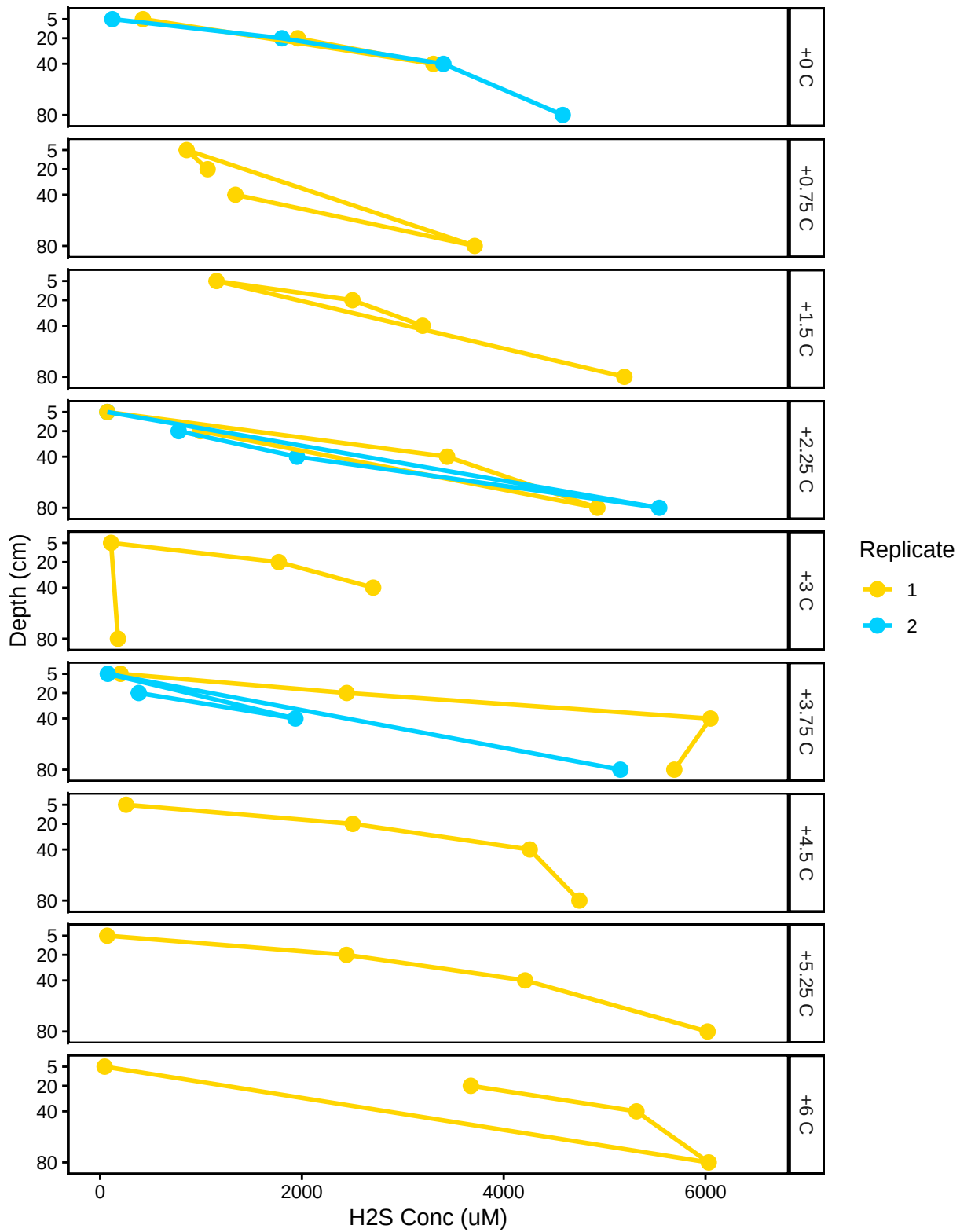

##Export Data

END